



A046

GYNECOLOGICAL TABLE

Registration code: SKF-B






SKSurgitech[™]
by Saikang Medical

Care you precisely

A046

FEATURES AND HIGHLIGHTS

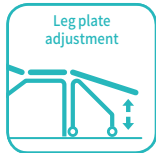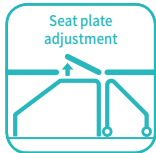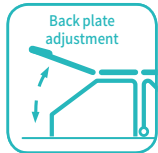

Technical parameters

▪ External size (LxWxH)	1900x600x800 mm
▪ Back plate up/down	7°/45°
▪ Seat plate up/down	8°/25°
▪ Leg plate up/down	25°/25°

Technical configuration

<input checked="" type="checkbox"/> Crank Wheel	1set
<input checked="" type="checkbox"/> Hand Bar	1set
<input checked="" type="checkbox"/> Anti-slip Bed Feet Cover	4pcs
<input checked="" type="checkbox"/> Mechanical Level	2sets
<input checked="" type="checkbox"/> Lithotomy Rod	1set
<input checked="" type="checkbox"/> Shoulder Holder	1set
<input checked="" type="checkbox"/> Stirrup	1set
<input checked="" type="checkbox"/> Arm Holder	1set
<input checked="" type="checkbox"/> 3" Castors with Break	4pcs



Structural design

The bed is composed of two parts: main bed and leg bed. The main bed table is divided into two parts: the back plate and the seat plate.



Application

This delivery bed is suitable for general operations such as maternal delivery, prenatal examination, abortion surgery, gynecological examination, etc.



Bed frame

The main bed and the auxiliary bed are separated separately, or they can be connected into one.



Leg bed table

The upper and lower folds of the leg bed table are adjusted by the ratchet rack mechanism to ensure the comfort of the patient's legs during the operation and the main bed table folds to meet the needs of the posture of the operation.



Mattress

Composed of 3 pieces, connected by a push button, with a sponge inside the artificial leather, leather fireproof, yellowing resistance, anti-mildew ball, twists and turns, abrasion resistance 100,000 times.




Fixed and movable casters

4 non-slip foot pads, 4 movable casters, and The foot brake locking device is attached to ensure the stability of the leg bed when it is used, and it is more convenient to push.



Anesthesia infusion stand

Use to hang medicine bottle during infusion, height can be adjusted.

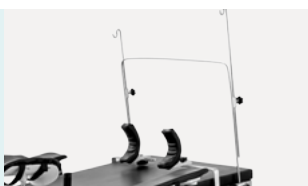

Shoulder Holder



Arm Holder



Lithotomy Rod



Stirrup



Hand Bar






CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@sksurgitech.com
🌐 www.sksurgitech.com



Company video



General catalogue