

# A100-4


## ELECTRIC OPERATING TABLE

Registration code: SKL-A



# A100-4

## FEATURES AND HIGHLIGHTS

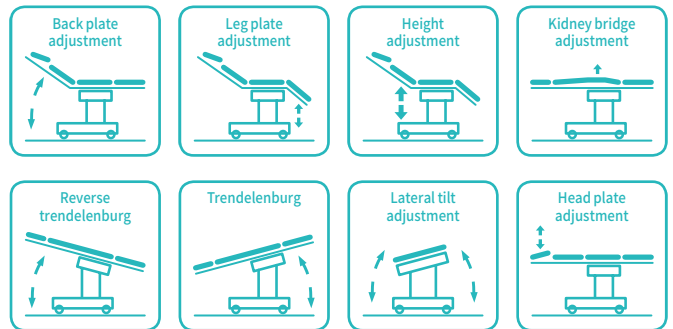

### Technical parameters

External Size (LxWxH)	2100x500x800-1000 mm
Built-in Kidney Bridge	0-110 mm
Back Plate Up/Down	75°/20°
Head Plate Up/Down	45°/90°
Leg Plate Up/Down	15°/90°
Trendelenburg/Reverse Trendelenburg	22°/25°
Lateral Tilt Left/Right	20°/20°
Horizontal Sliding	300 mm

### Technical configuration

<input checked="" type="checkbox"/> Bed Frame	1pc
<input checked="" type="checkbox"/> Anesthesia Screen Frame	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Kidney Bridge Handle	1pc
<input checked="" type="checkbox"/> Arm Support with Cushion	2pcs
<input checked="" type="checkbox"/> Leg Support with Cushion	2pcs
<input checked="" type="checkbox"/> Shoulder Support with Cushion	2pcs
<input checked="" type="checkbox"/> Body Support with Cushion	2pcs
<input checked="" type="checkbox"/> Foot Support with Cushion	2pcs

#### Remote Panel

All function controlled by hand controller, and can be locked.



#### C-arm

Extra-long horizontal sliding, available for c-arm operation.




#### Kidney Bridge

Built-in kidney bridge, more convenient and nice appearance.



#### Wash Basin Hole

Five section table top including Head, Back, Shoulder, Hip and Leg, all section applicable for X ray use and under each section with X ray cassette layer.

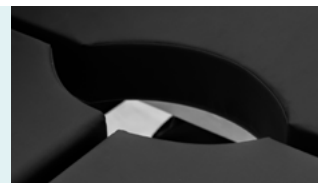

#### Head Plate

Detachable double head plate and leg plate meet operation position needs and use Germany cylinder assist system.



#### Powder Supply

Internal battery maintaining power supply for 50-80 operations and with low battery warning function.



Shoulder Support with Cushion



Body Support with Cushion



Arm Support with Cushion



Leg Support with Cushion



Lowest Position of Table Top  
**800mm**



Highest Position of Table Top  
**1000mm**



Longitudinal Adjustment  
**0-300mm**



Longitudinal Adjustment  
**0-300mm**



Reversed Trendelenburg  
**0°-22°**



Trendelenburg  
**0°-25°**



Lateral Tilt Left  
**0°-20°**



Lateral Tilt Right  
**0°-20°**



Kidney Bridge Elevation  
**110mm**



Anorectum Operation  
Position



E.N.T Operation Position



(Optional) Head Rest



(Optional) Lateral Arm  
Support



(Optional) Waste Basin For  
Gynecology



(Optional) Orthopedics Traction  
Frame





## **CONSTANT INNOVATION**



Design



Research



Manufacturing



Maintenance



Service

## **WORLDWIDE BRANCHES**

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa






**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
 +86 18021231901 ·  24/7  
 export@sksurgitech.com  
 www.sksurgitech.com



Company video



General catalogue