



A105

OBSTETRIC TABLE

Registration code: SKF-C




SKSurgitech[™]
by Saikang Medical

Care you precisely

A105

FEATURES AND HIGHLIGHTS


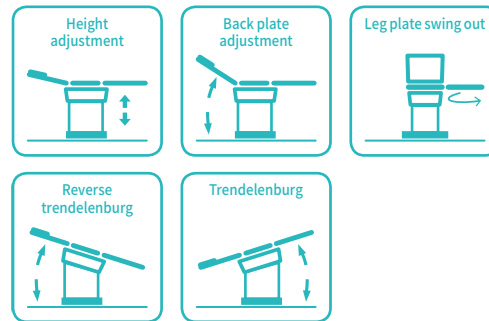


Technical parameters

External Size (LxWxH)	1900x600x700-900 mm
Back Plate Up/Down	50°/3°
Trendelenburg/Reverse Trendelenburg	10°/20°
Leg Board Outwards	90°
Working Area	600x600 mm
Patient Weight Capacity	185kg


Technical configuration

<input checked="" type="checkbox"/> Bed Frame	1pc
<input checked="" type="checkbox"/> Mattress	1pc
<input checked="" type="checkbox"/> Hand Bars	2pcs
<input checked="" type="checkbox"/> Arm Support with Cushion	2pcs
<input checked="" type="checkbox"/> Leg Support with Cushion	2pcs
<input checked="" type="checkbox"/> Anesthesia Screen Frame	1pc
<input checked="" type="checkbox"/> Stainless Steel Basin	1pc
<input checked="" type="checkbox"/> Handwheel	4pcs
<input checked="" type="checkbox"/> Oil Pot (With Oil)	1pc

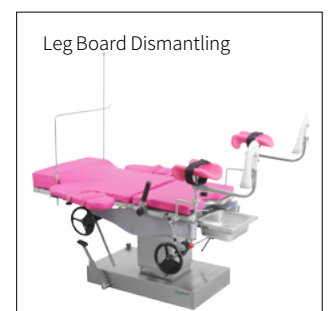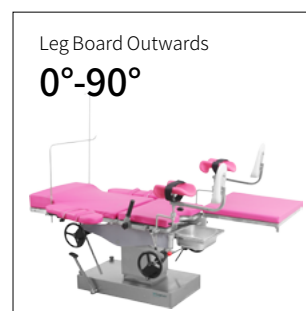

Material

The base seat and upright column covering are made of high quality 304 stainless steel, which is anti-corrosive and can be easily cleaned.
Optional: Powder Coated Steel Base.

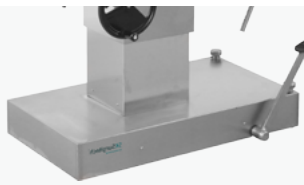

Function

The main bed and auxiliary bed are connected with each other, auxiliary table is easy for rotating and dismantling.



Rotary Brake Device

Rotary brake device which is easy for moving obstetric table.



Foot Step Hydraulic Elevating

Foot step hydraulic elevating, Imported Y type sealing ring and is durable.



Mattress

The mattress adopts high density memory sponge and is once formed, with no seam.
Optional: Detachable mattress.



Stainless Steel Sewage Basin

Equipped with stainless steel sewage basin and other various accessories.



Leg Support with Cushion

Double using leg support, can be used as boatstretcher.



Anesthesia Screen Frame



Arm Support with Cushion



Hand Bars



Handwheel





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690



Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

 +86 18021231901 ·  24/7

 export@sksurgitech.com

 www.sksurgitech.com



Company video



General catalogue