



# A106-2

## ELECTRIC-HYDRAULIC OPERATING TABLE

Registration code: SKL-A



**SKSurgitech**<sup>™</sup>  
by Saikang Medical

Care you precisely

# A106-2

## FEATURES AND HIGHLIGHTS


### Technical parameters

▪ External Length	2080mm
▪ External Width	550mm
▪ Table Height Min.	500mm
▪ Table Height Max.	800mm
▪ Reverse Trendelenburg	0°-25°
▪ Trendelenburg	0°-25°
▪ Lateral Tilt Left/Right	0°-20°/0°-20°
▪ Head Plate Up/Down	0°-45°/0°-90°

▪ Back Plate Up/Down	0°-80°/0°-20°
▪ Foot Plate Up/Down	0°-15°/0°-90°
▪ Foot Plate Horizontal Positioning	0°-90°
▪ Patient Weight Capacity	250kg

### Technical configuration

<input checked="" type="checkbox"/> Anaesthesia Screen	1pc
<input checked="" type="checkbox"/> Hand Remote	1pc
<input checked="" type="checkbox"/> Kidney Bridge Handle	1pc
<input checked="" type="checkbox"/> Arm Support	1pc
<input checked="" type="checkbox"/> Foot Support	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Shoulder Support	1pair
<input checked="" type="checkbox"/> Long Fixing Clamp	1pair
<input checked="" type="checkbox"/> Fixing Clamp	7pcs
<input checked="" type="checkbox"/> Power Line	1pc

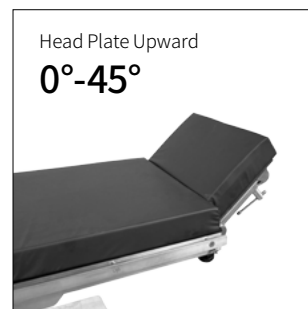

**Material**

The base, column and side hoder are all made of high quality 304 stainless steel, which is anti-corrosive and can be easily cleaned.



**Powder Supply**

Internal battery maintaining power supply for 50-80 operations and with low battery warning function.




**Remote Panel**

All function controlled by hand controller in the column, and can be locked.




**Hydraulic System**

Brand hydraulic system (German Hoerbiger) make table with low noise, high safety and reliability.

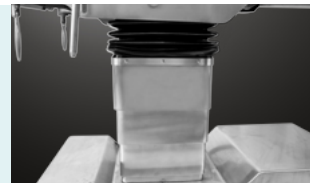


**Control System**

Foot treadle brake device which is easily for moving operating table.

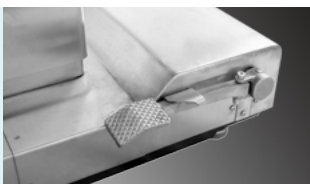

**Super Low Position**

Super low position (500mm) , suitable for operation on head.

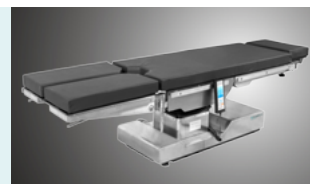

**Kidney Bridge**

Built-In Kidney Bridge, More Convenient And Nice Appearance.



**Screen**



**Armrest**




**Knee Curtch**



**Shoulder Support**



**Body Support**

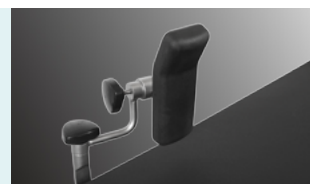




## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai,  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎️ +86 18021231901 · 🕒 24/7  
✉️ export@sksurgitech.com  
🌐 www.sksurgitech.com



Company video



General catalogue