



A2000

ELECTRIC OPERATING TABLE

Registration code: SKL-A



SKSurgitechTM
by Saikang Medical

Care you precisely

A2000

FEATURES AND HIGHLIGHTS



Technical parameters


External size (LxWxH):	2100x480x(680-930)mm
Reverse trendelenburg:	0-22°
Trendelenburg:	0-22°
Lateral tilt left/right:	0-15°
Head plate up/down:	55°/90°
Back plate up/down:	50°/15°
Leg plate up/down:	45°/90°
Leg plate outwards:	0-90°
Kidney bridge elevation:	0-100mm
Safe working load:	150kg
X-rays translucent:	NO
C-arm compatible:	NO
Multi-functional orthopedic frame:	YES

Technical configuration

<input checked="" type="checkbox"/> Bed frame	1 pc
<input checked="" type="checkbox"/> Mattress	1 set
<input checked="" type="checkbox"/> Anesthesia screen frame	1 pc
<input checked="" type="checkbox"/> Shoulder holder	2 pcs
<input checked="" type="checkbox"/> Arm holder	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Foot support	2 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Sliding block	10 pcs
<input checked="" type="checkbox"/> Spanner	4 pcs


Material

The base material for the table is high thickness 304 stainless steel, seamless welding by robot, with modern appearance and more stable.

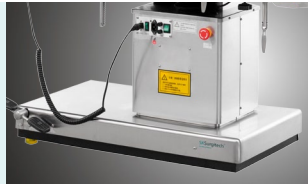

Mattress

Removable memory mattress, good rebound resilience, anti-static electricity and anti-microbial.

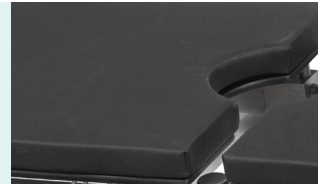

Brake

Pedal brake, easy to move and lock.



Safety

The table is provided with an emergency stop switch on the base column, help to stop the operating table running in case of emergency situation, ensure the safety of operation.



Leg plate

Leg plate with imported gas spring for easy control.




Kidney bridge

Manual kidney bridge, more convenience.



Battery could support

Battery could support about 50-80 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.

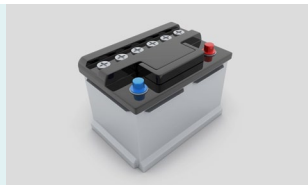

Armrest



Anesthesia Frame




Leg Crutches



Shoulder Support



Body Supporter

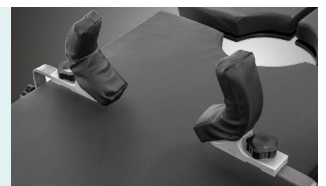

Professional functions

Lowest position of table top

680mm



Highest position of table top

930mm



Back plate up

0°-70°



Back plate down

0°-15°



Reversed trendelenburg

0°-25°



Trendelenburg

0°-25°



Right lateral tilt

0°-15°



Left lateral tilt

0°-15°



Head plate upward

0°-45°



Head plate downward

0°-90°



Leg section evagination

0°-90°



Leg plate up/down

-45°-90°



Kidney bridge elevation

0-100mm





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@sksurgitech.com
🌐 www.sksurgitech.com



Company video



General catalogue