



A202

ELECTRIC  
OPERATING TABLE

Registration code: SKL-A




SKSurgitech™  
by Saikang Medical

Care you precisely

# A202

## FEATURES AND HIGHLIGHTS

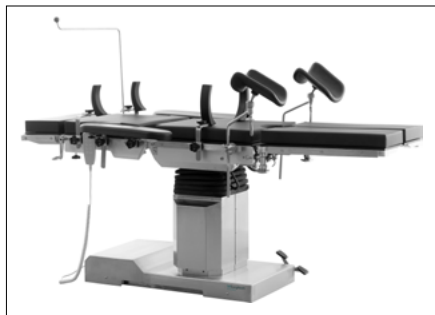

### Technical parameters


External size (LxWxH):	2000x500x(750-1050)mm
Reverse trendelenburg:	0-22°
Trendelenburg:	0-25°
Lateral tilt left/right:	0-20°
Head plate up/down:	45°/90°
Back plate up/down:	70°/35°
Leg plate down:	0-90°
Leg plate outwards:	0-90°
Kidney bridge elevation:	0-80mm
Horizontal sliding:	0-300mm
Safe working load:	360kg

### Technical configuration

<input checked="" type="checkbox"/> Arm support	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Body Support	2 pcs
<input checked="" type="checkbox"/> Shoulder Support	2 pcs
<input checked="" type="checkbox"/> Slider	7 pcs
<input checked="" type="checkbox"/> Tightener	2 pcs
<input checked="" type="checkbox"/> Anaesthetic Screen	1 pc
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Handwheel	1pc


**Material**

The chassis cover is made of high quality Ni-Cr stainless alloy steel by seamless robot welding, corrosion resistance, disinfection resistance, rust-proof and durable.

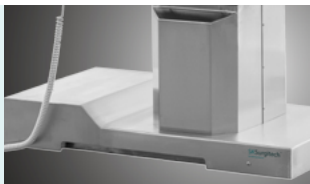

**Mattress**

Removable memory mattress, good rebound resilience, anti-static electricity and anti-microbial.

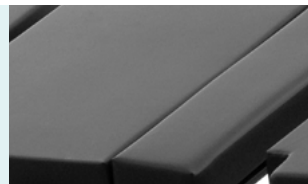

**Control System**

Electric-hydraulic control system, non-gear structure, all functions are controlled by the hand controller.



**Braking system**

Mechanical brake device at the side of the base, with fixed or moving function, flexible and reliable.



**Hand controller**

All functions are controlled by the hand controller.



**Kidney bridge**

Kidney bridge is adjusted by the hand wheel at the side of head board, more convenient, with separate mattress.



**Battery could support**

Battery could support about 50-80 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.



**Horizontal movement**

300mm horizontal movement function, with high strength phenolic resin table top, could be used for X-ray and C-arm.



**Leg plate**

The up and down folding of the leg plate is manually controlled by the high strength gear, which makes the operation more effort-saving and smooth.



**Armrest**



**Anesthesia Frame**



**Leg Crutches**



**Shoulder Support**



**Body Supporter**



Professional functions

Lowest position of table top  
**750mm**



Highest position of table top  
**1050mm**



Back plate up  
**0°-70°**



Back plate down  
**0°-35°**



Reversed trendelenburg  
**0°-22°**



Trendelenburg  
**0°-25°**




Right lateral tilt  
**0°-20°**



Left lateral tilt  
**0°-20°**



Head plate upward  
**0°-45°**

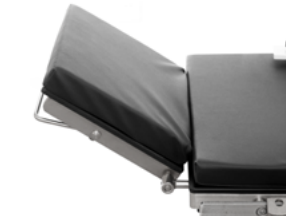

Head plate downward  
**0°-90°**



Leg section evagination  
**0°-90°**



Leg plate up/down  
**-45°-90°**



Kidney bridge elevation  
**0-80mm**



Horizontal movement  
**0-300mm**





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎️ +86 18021231901 · 🕒 24/7  
✉️ export@sksurgitech.com  
🌐 www.sksurgitech.com



Company video



General catalogue