



# A3008A-1

## MANUAL OPERATING TABLE

Registration code: SKF-C




**SKSurgitech**<sup>™</sup>  
by Saikang Medical

Care you precisely

# A3008A-1

## FEATURES AND HIGHLIGHTS

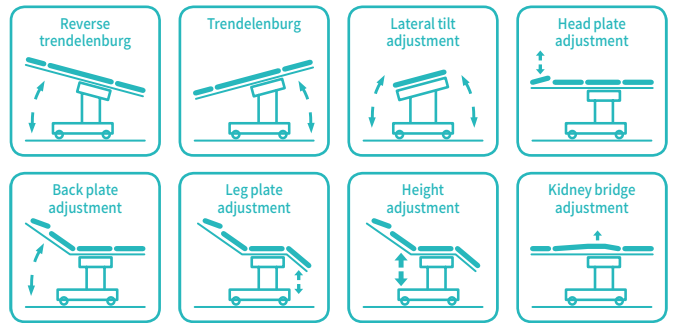

### Technical parameters

External Size (LxWxH):	2100x480x750-950 mm
Built-in Kidney Bridge:	0-100 mm
Back Plate Up/Down:	75°/5°
Head Plate Up/Down:	30°/90°
Trendelenburg/Reverse Trendelenburg:	30°/30°
Lateral Tilt Left/Right:	20°/20°
Leg Plate Down:	90°
Patient Weight Capacity:	250 kg
X-rays translucent:	optional with extra charge
C-arm compatible:	NO
Multi-functional orthopedic frame:	YES

### Technical configuration

<input checked="" type="checkbox"/> Anesthesia Screen Frame	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Kidney Bridge Handle	1pc
<input checked="" type="checkbox"/> Arm Support with Cushion	2pcs
<input checked="" type="checkbox"/> Leg Support with Cushion	2pcs
<input checked="" type="checkbox"/> Shoulder Support with Cushion	2pcs
<input checked="" type="checkbox"/> Body Support with Cushion	2pcs
<input checked="" type="checkbox"/> Long Fixing Clamp	1pair
<input checked="" type="checkbox"/> Fixing Clamp	8pcs



**Application**

This operating table is a comprehensive multi-function type, suitable for medical units to perform various operations, all movements are controlled by the manual cranks.



**Material**

The base, column and side holder are all made of high quality 304 stainless steel, which is anti-corrosive and can be easily cleaned.



**Control System**

Foot treadle brake device which is easily for moving operating table.



**Built-In Kidney Bridge**

Built-In kidney bridge, more convenient and nice appearance.



**Hydraulic**

The height adjustment controlled by hydraulic pump.




**Mattress**

PU cover, custom mattress, easy disassembly and cleaning, thickness of mattress: 6cm



**Shoulder Support with Cushion**

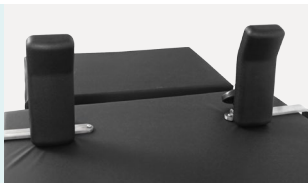

**Anesthesia Screen Frame**



**Arm Support with Cushion**



**Leg Support with Cushion**



**Body Support with Cushion**





## **CONSTANT INNOVATION**



Design



Research



Manufacturing



Maintenance



Service

## **WORLDWIDE BRANCHES**

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa





**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
 +86 18021231901 ·  24/7  
 export@sksurgitech.com  
 www.sksurgitech.com



Company video



General catalogue