



A98-3

ELECTRIC OBSTETRIC BED

Registration code: SKL-A





SAIKANG / Medical

Healthier
your life

A98-3

FEATURES AND HIGHLIGHTS

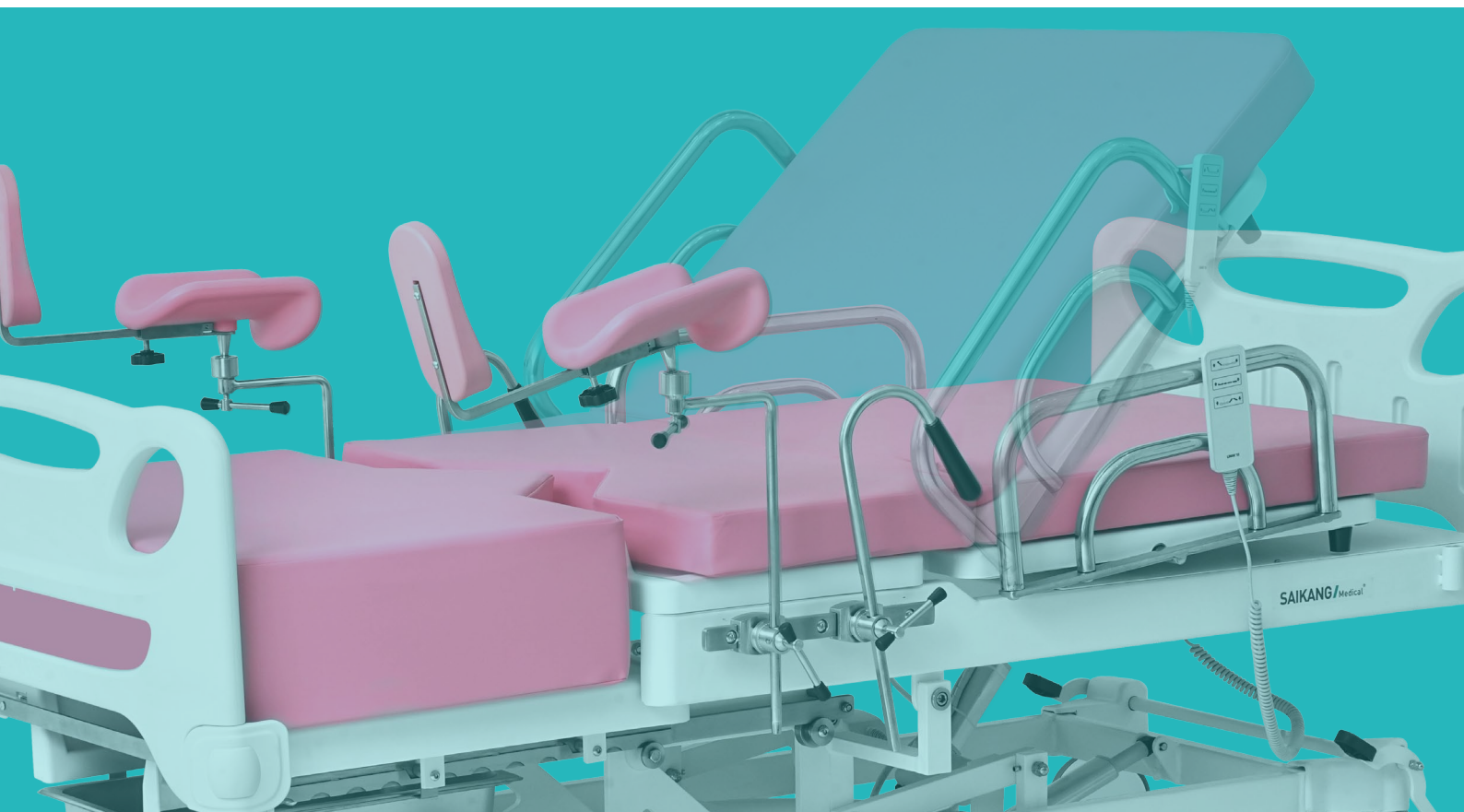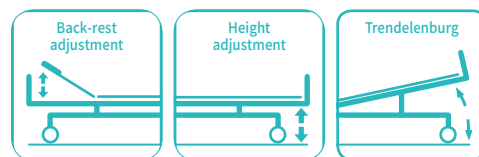

Technical parameters

▪ External size (LxW)	2150 x 980 mm
▪ Height Adjustment	590-915 mm
▪ Back-rest adjustment	0°-70°
▪ Trendelenburg	0°-10°

Technical configuration

<input checked="" type="checkbox"/> Pull Handle	2pcs
<input checked="" type="checkbox"/> Leg support	2pcs
<input checked="" type="checkbox"/> Foot support	2pcs
<input checked="" type="checkbox"/> Hand Remote Controller	1pc
<input checked="" type="checkbox"/> Stainless Steel Basin	1pc
<input checked="" type="checkbox"/> Instruction Manual	1pc
<input checked="" type="checkbox"/> 6" Castor	1pc
<input checked="" type="checkbox"/> Mattress	1set

Material

The base seat and height adjustment system are made of powder coated steel, high gloss, anti corrosion, easy to clean.

Bed board

Bed surface has PP cover, made by virgin material. The PP cover is detachable, and easy to clean.

Head&Foot Board

The bed platform and head&foot board made by virgin PP material, the head&foot board is movable.

Hand controller

Equipped with one hand controller.

Guardrail

Handrail is rotatable. It is convenient for the transfer of puerpera.

Handle

The hand bar can be dismantled and movable.

Lithotomy rod

The leg support can be dismantled and movable.

Foot support

The foot support can be dismantled and movable.

Castor

6 inch central controlled castor for moving obstetric table. TPR tire, quiet and wear-resistant, can be washed and sterilized by high temperature, save anti-winding hard shell, united forming without bolts.

Accessory

Equipped with stainless steel sewage basin and other various accessories.

Mattress

The mattress adopts high density memory sponge and is once formed. Mattress is movable.

Foot section

Foot section can be retracted into the bed, convenient for gynecological and delivery surgical.



External Size (LxW)

2150x980mm



Back-rest adjustment

70°



Height adjustment

590-915mm



Trendelenburg


10°



Foot section into the bed



CONSTANT INNOVATION

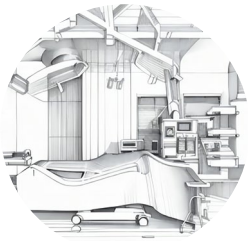


Design



Research



Manufacturing

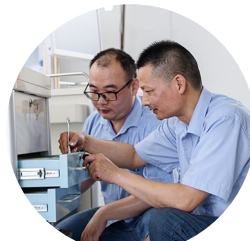

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited

Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC

17800 CASTLETON ST STE 180,
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue