



# A98-5

## ELECTRIC OBSTETRIC BED

Registration code: SKL-A




**SKSurgitech**<sup>™</sup>  
by Saikang Medical

Care you precisely

# A98-5

## FEATURES AND HIGHLIGHTS

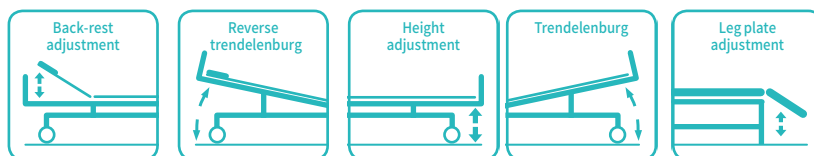

### Technical parameters

External size (LxW):	2140x1080x(750-900)mm
Height(without mattress):	580-790mm(±10mm)
Back-rest Adjustment:	0-70°
Reverse Trendelenburg:	0-6°
Trendelenburg:	0-6.5°
Overall height adjustment:	750mm-900mm ±10mm
Leg-rest Adjustment:	0-60°
Power supply:	AC220V
Input power:	250VA

### Technical configuration

<input checked="" type="checkbox"/> Mattress	1 pc
<input checked="" type="checkbox"/> Basin	1 pc
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Handle	2pcs
<input checked="" type="checkbox"/> Leg Support	2pcs
<input checked="" type="checkbox"/> Castors	4pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2pcs

# QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



## Design

The overall design is humanized and modern, and it is an obstetric and disease comprehensive bed integrating obstetric bed, hospital bed, examination bed and operating bed.



## Bed material

The bed frame is welded by 30\*60mm,  $\geq 1.5$ mm high-quality rectangular steel, beautiful and burr-free, with strong bearing capacity. It can bear up to 200kg.



## Paintwork

The surface of the bed body is treated with electrostatic spraying, full anti-rust treatment, anti-aging, not easy to rust, antibacterial and environmentally friendly.



## Welding | Craftsmanship

The overall welding is adopted by the robot, and the appearance is smooth, stable and reliable.



## Back-rest

The back section can be raised and lowered to assist pregnant women to sit and stand, and relieve the pressure of pregnant women during labor.



## Leg-rest

The leg section can be lifted and lowered to assist pregnant women to give birth smoothly, and the posture can be adjusted arbitrarily, such as sitting, kneeling, lying down, and squatting.



## Trendelenburg&Reverse Trendelenburg

Electric control is used for trendelenburg and reverse trendelenburg.



## Height Adjustment

The overall height of the bed body can be adjusted to facilitate the adjustment of the height of the bed, which is more suitable for nursing staff to take care of pregnant women.

**Guardrail**

Imported PP one-time molding guardrail, double hydraulic rods, beautiful and easy to clean.

**Head & foot board**

European-style head and foot board, resistant to fading and aging, high strength, easy to clean, firm and strong; hidden locking switch, stable and reliable, easy to disassemble, more suitable for the concept of a family-oriented and warm delivery room.

**Motor and controller**

Brand L&K DC motor lifting system, constant speed, quiet, no static, safe voltage (24V).

**Caster**

It adopts high-grade silent polyurethane wear-resistant casters with a diameter of  $125 \pm 5$ mm, which are silent and wear-resistant, and move smoothly. The built-in fully enclosed self-sliding bearing; one foot brake, four wheel brake.

**Mattress**

High-density foam mattress, not removable. The U-shaped incision design at the buttocks is more convenient for doctors to operate.

**Leg support**

Manual multi-angle adjustment, more convenient to use.

**Bed surface**

Made of high-grade leather, it is comfortable and beautiful. The pink bed surface design can relieve the psychological pressure of parturients; the 8cm thick high-density sponge is comfortable and elastic.

**Handle**

The power-assisted handle can be adjusted at multiple angles, and the humanized design is safer and more convenient.

**Central locking system**

The side central control brake system makes the delivery bed more stable.

**Basin**

The stainless steel dirt basin is standard, strong and durable, easy to clean, easy to install and disassemble.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi,  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai,  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎️ +86 18021231901 · 🌐 24/7  
✉️ export@sksurgitech.com  
🌐 www.sksurgitech.com



Company video



General catalogue