



A98-8

ELECTRIC OBSTETRIC BED

Registration code: SKL-A



SAIKANG Medical

Healthier
your life

A98-8

FEATURES AND HIGHLIGHTS



Technical parameters

▪ Bed Dimensions (LxW)	2350 x 1080 mm
▪ Height Adjustment	470-880 mm
▪ Bed Platform Dimensions (LxW)	1980 x 900 mm
▪ Back-rest adjustment	0°-65°
▪ Seat plate lifting	15°
▪ Trendelenburg	8°
▪ Foot plate height adjustment	0-200mm
▪ Footboard Swing Out	0-90°
▪ Footboard Uplift	0-90°
▪ Mattress Thickness	0-130mm
▪ Power Supply	AC220V
▪ Power Input	250VA
▪ Packing Dimensions	1800*950*880mm
▪ Weight	N.W.: 200KGS, G.W.: 270KGS

Technical configuration

<input checked="" type="checkbox"/> Guardrail	2pcs
<input checked="" type="checkbox"/> Head&foot board	1pc
<input checked="" type="checkbox"/> Motor	1set
<input checked="" type="checkbox"/> Obstetric handle	2pcs
<input checked="" type="checkbox"/> Auxiliary table	1pc
<input checked="" type="checkbox"/> Delivery holder	1pc
<input checked="" type="checkbox"/> Foot support	1pc
<input checked="" type="checkbox"/> 6" Castor, Central locking system	4pcs
<input checked="" type="checkbox"/> Wall bumper wheel	2pcs
<input checked="" type="checkbox"/> Night light	1pc
<input checked="" type="checkbox"/> Mattress	1set
<input checked="" type="checkbox"/> Basin	1pc
<input checked="" type="checkbox"/> IV pole	1pc
<input checked="" type="checkbox"/> Backup battery	1pc

Function

- Super low starting position facilitates easy access for patients to get on and off the bed.
- Backrest elevation Assists pregnant women in sitting and standing, alleviating the stress of childbirth.
- Folding hip plate is convenient for childbirth.
- Manual Trendelenburg Operation is easy to use for adjusting the bed's position.
- Automatic seat tilting function allows for easy operation.
- The bed is designed with rebound protection for the lowest and trendelenburg position.
- Equipped with an audio system and USB port, allowing medical staff and parturients to select their preferences for stress relief.
- One-button lock to prevents accidental operation and effectively protects the safety of the parturient.

Design

It is a humanised maternity bed under the guidance of the new idea of international accompanied delivery, which integrates the reception of labour, waiting for labour, delivery, recuperation and general obstetric surgery, general examination and diagnosis. The overall design is humanised and can be used for delivery in supine position, sitting position, side-lying position, squatting position and other positions.



Bed Height Min.

470 mm



Bed Height Max.

470 mm



Back-rest Adjustment

65°



Seat Plate Lifting

10°



Foot Plate Height Adjustment

0-200 mm



Footboard Swing Out/Footboard Uplift

0-90°/0-90°

Material

The whole bed frame is a steel-plastic structure, the maximum thickness of steel can be up to 10mm
Load capacity: 270KG.

Guardrail

Split lifting design guardrail, PP material (or the same level of material) one-time molding, the surface is soft to provide caregivers and maternal health with warm and comfortable feeling. Nurse and maternal control panels are embedded on both sides, convenient for both sides to adjust different maternal positions.

Head&Foot Board

The removable headboard in PE plastic or equivalent material fits the concept of a caring and cosy maternity ward.

Bed platform

The bed platform is blow-moulding by moulding tools, with anti-corrosion, anti-odour, easy to clean.
Load capacity ≥ 150 KG.

Surface coating

Bed surface electrostatic spraying treatment, using imported plastic powder, durable, anti-rust, antibacterial and environmental protection.

Motor

Timosun linear motor, no need for regulator, safe and smooth operation without noise. DC motor lifting system, constant speed, silent, no static electricity
safe voltage : DC24V.

Obstetric handle

Folded handle with lock, ergonomic design, safer and more convenient. The handles can be moved back and forth to adjust 150mm, which is convenient for patients.

Auxiliary table and delivery holder

Standardly equipped with large and small auxiliary table and delivery holder for free positions, which is convenient for patients to choose different positions for delivery.

Emergency CPR button


Emergency CPR button, to improve the speed of emergency resuscitation.

Foot support

Foot support with one-handed operation to adjust the maternal feet to the appropriate position; leg rest frame is raised or lowered with legrest lifting button.

Central locking brake system

Central locking brake system makes the stability of the delivery bed better, with directional and universal function.

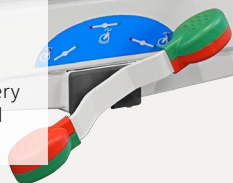

Castor

6-inch high-quality medical castors, double-wheel design, to prevent sticking in of hair.




Wall bumper wheel

The bed is equipped with wall bumper wheel to effectively protect the bed and the equipment in the room.



Night light

Sensor-type night light provides automatic illumination of the lower bed area and the ground, and the light intensity can be adjusted in steps.

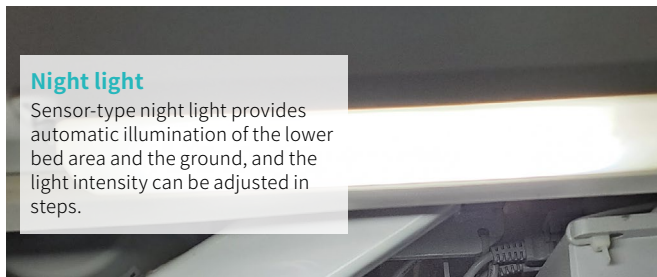

Mattress

Made of memory foam, the mattress is comfortable, soft, anti-corrosive, odor-resistant, and easy to clean. It ensures proper support for the parturient's waist and hips while alleviating pain and increasing comfort.




Basin

ABS Splash-Proof Basin, durable and stain-resistant, easy to install and remove.



IV pole

The main body of IV pole is made of stainless steel, and the height is adjustable to ensure the safety and reliability.

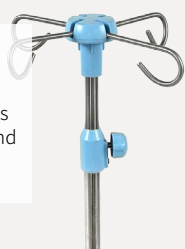


Backup battery

The uninterrupted power supply system ensures use during emergencies and features a display for battery life and voltage (capable of continuous use for 48 hours).





CONSTANT INNOVATION

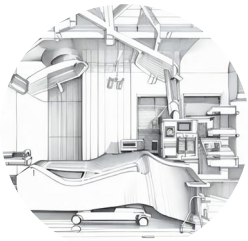


Design



Research



Manufacturing

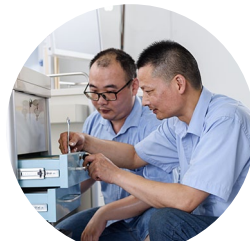

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
17800 CASTLETON ST STE 180,
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue