



# A99-12

## ELECTRIC OBSTETRIC BED

Registration code: SKF-C



**SKSurgitech**<sup>™</sup>  
by Saikang Medical

Care you precisely

# A99-12

## FEATURES AND HIGHLIGHTS



### Technical parameters

▪ External (LxW):	1800x600mm
▪ Table height:	700-950mm
▪ Back plate up folding:	0-60°
▪ Reverse Trendelenburg:	0-25°
▪ Trendelenburg:	0-16°
▪ Load bearing:	≤150kg




### Technical configuration

<input checked="" type="checkbox"/> Arm support	2 pcs
<input checked="" type="checkbox"/> Handle	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Stirrup	2 pcs
<input checked="" type="checkbox"/> Slider	4 pcs
<input checked="" type="checkbox"/> Tightener	2 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Foot switch	1 pc
<input checked="" type="checkbox"/> Stainless steel basin	1 pc



**Material**

The chassis cover is made of high quality 304 stainless steel, which is welded by the robot seamlessly. The shape is beautiful and stable.

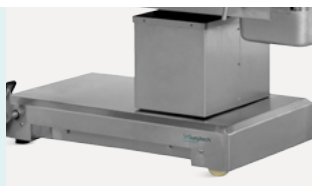

**Mattress**

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.



**Hand controller**

All the functions are controlled by hand controller.



**Braking system**

The pedal controls the brake function of the casters.



**Control System**

Electric control the height lifting, trendelenburg and reverse trendelenburg, back plate folding up, easy and convenient operation.




**Battery**

The delivery bed is equipped with a standard battery, which could support the electric operation for 50-80 times, and could keep maintaining the operation of various functions required by the operator under special circumstances such as power failure.



**Pull handles**

Both sides of the delivery bed are equipped with collapsible handle, which is convenient for the puerpera to get on and off the bed, with anti-slip lines, to prevent slipping when the force is applied.

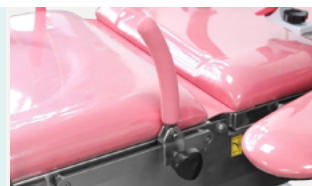

**Sewage basin**

The delivery bed is equipped with a stainless steel basin, which can be stored under the bed through a high-quality three-stage slide rails.



**Auxiliary Table**

The auxiliary table can be retractable and stored under the bed body, easy to use.



**Leg supports**

The leg support is integrated with stirrup function, and the delivery bed is equipped with a pair of stirrups.



**Arm Support**

The arm support could be adjusted into multiple angles, it can be used to comfortably fix the puerpera's arm at various angles.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
📞 +86 18021231901 · 🕒 24/7  
✉️ export@sksurgitech.com  
🌐 www.sksurgitech.com



Company video



General catalogue