



# A99-13

## ELECTRIC OBSTETRIC BED

Registration code: SKF-C




**SKSurgitech**<sup>™</sup>  
by Saikang Medical

Care you precisely


# A99-13

## FEATURES AND HIGHLIGHTS

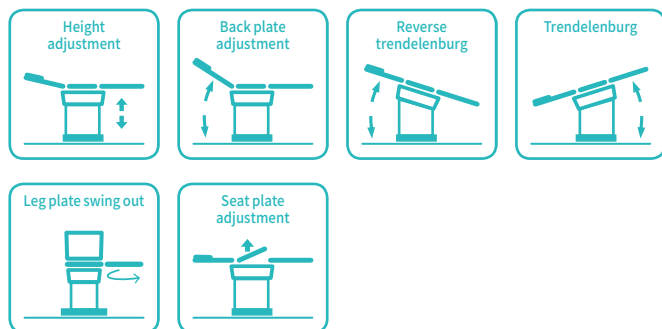

## Technical parameters

▪ External (LxW):	2100x600mm
▪ Table height:	750-1000mm
▪ Trendelenburg:	0-22°
▪ Reverse trendelenburg:	0-22°
▪ Back plate up:	0-60°
▪ Seat plate up:	0-16°
▪ Leg board outwards:	0-90°
▪ Arm board outwards:	0-180°
▪ Load bearing:	≤300kg

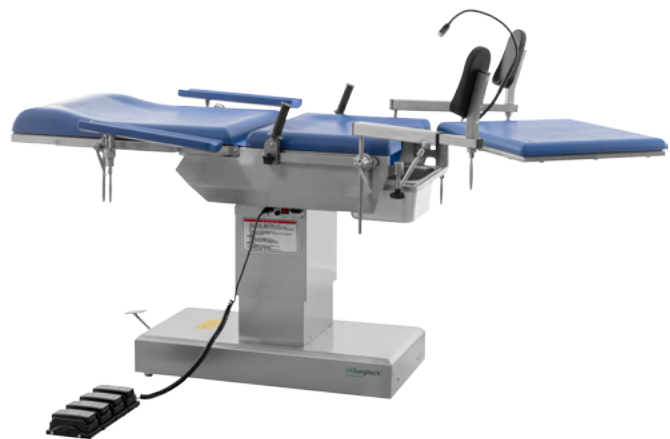


## Technical configuration

<input checked="" type="checkbox"/> Armrest	2 pcs
<input checked="" type="checkbox"/> Knee Curtch	2 pcs
<input checked="" type="checkbox"/> Pull handle	2 pcs
<input checked="" type="checkbox"/> Anesthesia frame	1 pc
<input checked="" type="checkbox"/> SST rinsing basin	1 pc
<input checked="" type="checkbox"/> Pedal controller	1 pc
<input checked="" type="checkbox"/> LED examination light	1 pc



**Material**

Base structure is made of high-quality 304 stainless steel, robot welded with high strength and modern appearance.

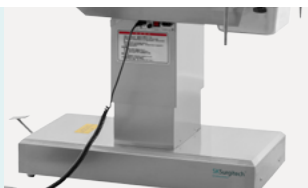

**Mattress**

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.



**Foot pedal controller**

All movements are controlled by foot pedal controller.



**Brake**

Pedal brake, easy to move and lock. Pedal down to release the inner castors down. Pedal up to retract the castors, which enables the table locked on the floor.



**LED examination light**

Standard accessory: LED examination light. More convenient to do gynecological examination.



**Arm plate**

Flexible arm plate adjustment.



**Battery could support**

Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table runing in the event of power failure.



**Leg supports**

Equipped with two different leg supports to meet different surgical needs.



**Auxiliary plate**

Auxiliary plate is easy for rotating and dismantling. Rotating angle: 90°.



**Pull handles**

Pull handles can be folded to horizontal position. Easy patient mobilisation.



**Sewage basin**

Equipped with stainless steel sewage basin and other various accessories.



**Electrical controlled**

All movements are electrical controlled: Back plate, seat plate, table elevation, trendelenburg and reverse-trendelenburg.






## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SKSurgitech™

by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
☎️ +86 18021231901 · 🕒 24/7  
✉️ export@sksurgitech.com  
🌐 www.sksurgitech.com



Company video



General catalogue