



# A99-5

## ELECTRIC OBSTETRIC BED

Registration code: SKF-C

**SKSurgitech™**  
by Saikang Medical

Care you precisely

# A99-5


## FEATURES AND HIGHLIGHTS

### Technical parameters

|                            |            |
|----------------------------|------------|
| ▪ External (LxW):          | 1950x660mm |
| ▪ Table height:            | 680-940mm  |
| ▪ Seat plate up folding:   | 0-20°      |
| ▪ Seat plate down folding: | 0-9°       |
| ▪ Back plate up folding:   | 0-45°      |
| ▪ Back plate down folding: | 0-10°      |
| ▪ Leg plate abduction:     | 0-90°      |
| ▪ Load bearing:            | ≤300kg     |

### Technical configuration

|   |       |
|---|-------|
| <input checked="" type="checkbox"/> Arm support           | 2 pcs |
| <input checked="" type="checkbox"/> Handle                | 2 pcs |
| <input checked="" type="checkbox"/> Leg support           | 2 pcs |
| <input checked="" type="checkbox"/> Stirrup               | 2 pcs |
| <input checked="" type="checkbox"/> Slider                | 4 pcs |
| <input checked="" type="checkbox"/> Tightener             | 2 pcs |
| <input checked="" type="checkbox"/> Hand controller       | 1 pc  |
| <input checked="" type="checkbox"/> Foot switch           | 1 pc  |
| <input checked="" type="checkbox"/> Stainless steel basin | 1 pc  |

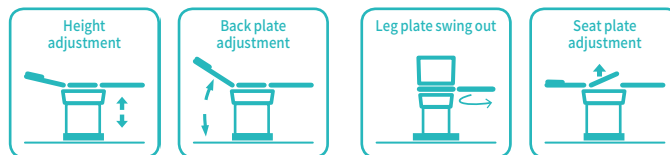

**Material**

The chassis cover is made of high quality 304 stainless steel, which is welded by the robot seamlessly. The shape is beautiful and stable.



**Mattress**

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.



**Foot pedal controller**

All the functions are controlled by hand controller and foot switch.



**Braking system**

The pedal controls the brake function of the casters.



**Control System**

Electric control the height lifting, seat plate up and down, back plate up and down, convenient operations.



**Arm Support**

The arm support could be adjusted into multiple angles with the straps, it can be used to comfortably fix the puerpera's arm at various angles.



**Leg supports**

The leg support is integrated with stirrup function, and the delivery bed is equipped with a pair of stirrups.




**Auxiliary Table**

The auxiliary table can be folded 90° outwards, used as a instrument table, detachable, and adapt to a variety of surgical environments.



**Pull handles**

Pull handles can be folded to horizontal position. Easy patient mobilisation.

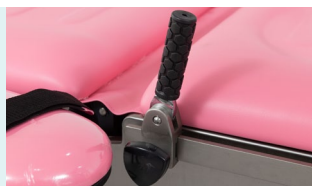

**Sewage basin**

Equipped with stainless steel sewage basin and other various accessories.





## **CONSTANT INNOVATION**



Design



Research



Manufacturing



Maintenance



Service

## **WORLDWIDE BRANCHES**

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa





**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
 +86 18021231901 ·  24/7  
 export@sksurgitech.com  
 www.sksurgitech.com



Company video



General catalogue