



A99-5

ELECTRIC OBSTETRIC BED

Registration code: SKF-C

SKSurgitech[™]
by Saikang Medical

Care you precisely

A99-5


FEATURES AND HIGHLIGHTS

Technical parameters

▪ External (LxW):	1880x600mm
▪ Table height:	680-940mm
▪ Seat plate up folding:	0-20°
▪ Seat plate down folding:	0-9°
▪ Back plate up folding:	0-45°
▪ Back plate down folding:	0-10°
▪ Leg plate abduction:	0-90°
▪ Load bearing:	≤300kg

Technical configuration

<input checked="" type="checkbox"/> Arm support	2 pcs
<input checked="" type="checkbox"/> Handle	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Stirrup	2 pcs
<input checked="" type="checkbox"/> Slider	4 pcs
<input checked="" type="checkbox"/> Tightener	2 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Foot switch	1 pc
<input checked="" type="checkbox"/> Stainless steel basin	1 pc

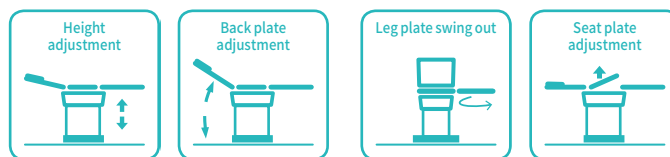

Material

The chassis cover is made of high quality 304 stainless steel, which is welded by the robot seamlessly. The shape is beautiful and stable.




Mattress

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.



Foot pedal controller

All the functions are controlled by hand controller and foot switch.

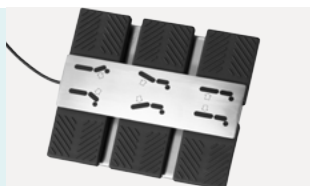

Braking system

The pedal controls the brake function of the casters.



Control System

Electric control the height lifting, seat plate up and down, back plate up and down, convenient operations.



Arm Support

The arm support could be adjusted into multiple angles with the straps, it can be used to comfortably fix the puerpera's arm at various angles.



Leg supports

The leg support is integrated with stirrup function, and the delivery bed is equipped with a pair of stirrups.



Auxiliary Table

The auxiliary table can be folded 90° outwards, used as a instrument table, detachable, and adapt to a variety of surgical environments.



Pull handles

Pull handles can be folded to horizontal position. Easy patient mobilisation.



Sewage basin

Equipped with stainless steel sewage basin and other various accessories.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@sksurgitech.com
🌐 www.sksurgitech.com



Company video



General catalogue