



A99-7

ELECTRIC OBSTETRIC BED

Registration code: SKL-A





SAIKANG / Medical

Healthier
your life

A99-7

FEATURES AND HIGHLIGHTS

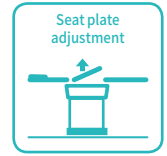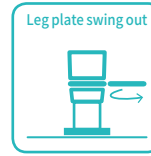

Technical parameters

▪ External (LxW):	1980x720mm
▪ Table height:	630-880mm
▪ Back plate up/down:	0-55°/0-15°
▪ Seat plate up/down:	0-35°/0-3°
▪ Leg board outwards:	0-90°

Technical configuration

<input checked="" type="checkbox"/> Bed Frame	1 pc
<input checked="" type="checkbox"/> Foot Pedal	1 pc
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Anesthesia Screen	1 pc
<input checked="" type="checkbox"/> Hand Bar	2 pcs
<input checked="" type="checkbox"/> Leg Support with Cushion	2 pcs
<input checked="" type="checkbox"/> Arm Support with Cushion	2 pcs
<input checked="" type="checkbox"/> Stainless Steel Sewage Basin	1 pc
<input checked="" type="checkbox"/> Castors	4 pcs

Accessories

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



Pull Handle




Armrest




Knee Curutch



Pedal Remote Controller

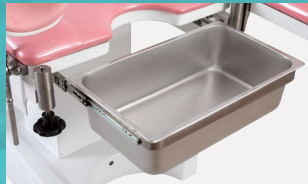

Stainless Steel Basin

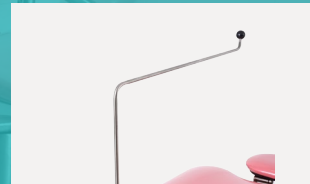

Anesthesia Frame



Height adjustment

630-880_{mm}



Back-rest up

55°



Back-rest down

18°



Seat plate up/down

0-35°/0-5°




Reverse trendelenburg

0-90°

QUALITY CONTROL

Material

The base seat is plastic cover and height adjustment system is made of iron, environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

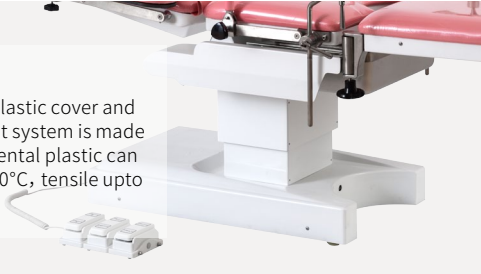


Functions

The main bed and auxiliary bed are connected with each other, auxiliary table is easy for rotating and dismantling.




Handle controller

Handle controller which is easy for moving obstetric table.

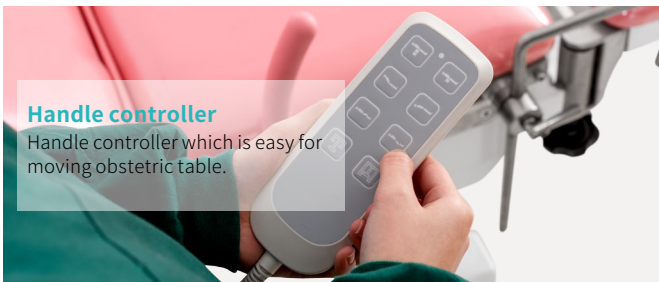


Leg support

Double using leg support, can be used as boatstretcher.

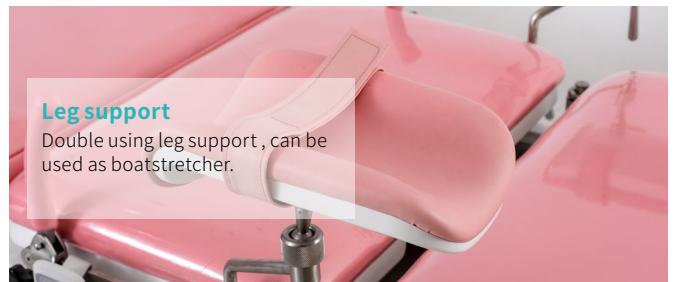


Mattress

The mattress adopts high density memory sponge and is once formed, with no seam.

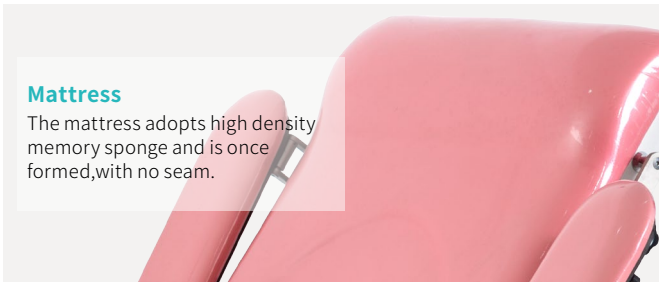

Bed width

seamless foam mattress, with waterproof PU cover. Easy cleaning and disinfection.

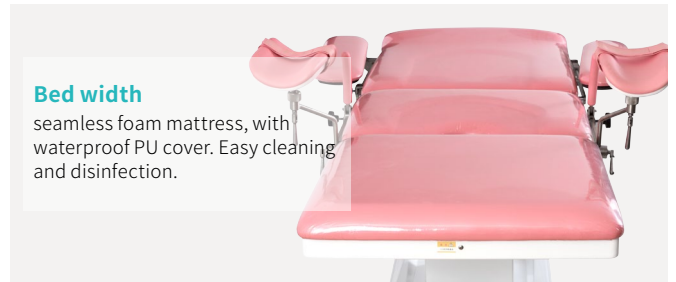


Control system

Good quality domestic motor, Various actions are controlled by remote controller.




Sewage basin

Equipped with stainless steel sewage basin and other various accessories.

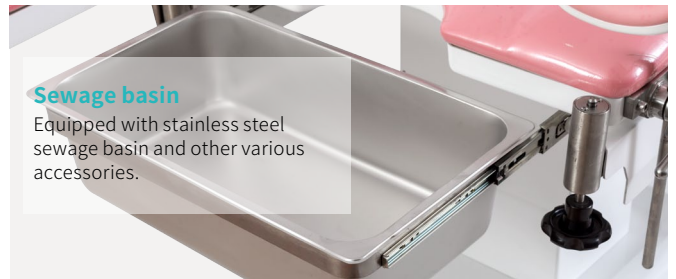



CONSTANT INNOVATION

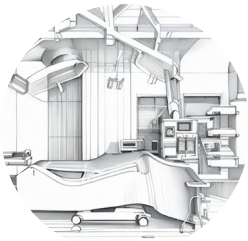


Design



Research



Manufacturing

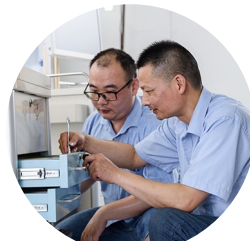

Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
 W3400 Corrientes, Argentina
 TEL: +55 (11) 2758-0438
 Rua Universal, Hollywood,
 São Bernardo do Campo,
 Brazil
 TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
 Mezzanine Floor, Taj Tower,
 upper hill, nairobi, Kenya
 TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
 Sala Thammasop Subdistrict
 Thawi
 Watthana District
 Bangkok Thailand
 TEL: +66 642741046

UK Branch

Uk Healthward International Limited
 DEPT 111,
 196 High Road Wood Green,
 London, N22 8HH
 United Kingdom
 TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
 454004, Russian Federation,
 Chelyabinsk,
 Universitetskaya
 Naberezhnaya, Russia
 TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
 Musaffah
 P.O Box 53200, Abu Dhabi,
 UAE
 Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
 STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue