



SK-L032

EXAMINING LAMP



SAIKANG / Medical

Healthier
your life



SK-L032

FEATURES AND HIGHLIGHTS

Suitable for clinical operation lighting



Technical parameters

▪ Diameter of lighthouse(mm):	250mm±10
▪ Illumination(at 1m)(Lux):	≥40000
▪ Brightness control range:	4levels adjustment
▪ Color temperature (Kelvin):	4000-5000
▪ Adjustable color temperature:	2levels adjustment
▪ Color rendering index(CRI):	≥90
▪ Bulb type:	LED
▪ Bulb brand:	OSRAM
▪ Bulb quantity:	7
▪ Max power(w):	12
▪ Average bulb life(hours):	50000
▪ Overall lamp height:	1.25-1.75m adjustment
▪ Power cable length:	2.5m
▪ Power supply voltage:	AC100-240V, 50/60HZ

Technical configuration


<input checked="" type="checkbox"/> 2-inch TPU universal silent wheels	4 pcs
<input checked="" type="checkbox"/> Base	1 set
<input checked="" type="checkbox"/> Control panel	1 set

Options

<input type="checkbox"/> Back-up battery	
--	--

Material technology

The lampshade shell is made of aluminum + integrated ABS structure design, which facilitates heat dissipation. The surface adopts the high-voltage electrostatic spraying process, and the imported environmentally friendly antibacterial plastic powder that used to ensures that the product meets the requirements of surgical hygiene.

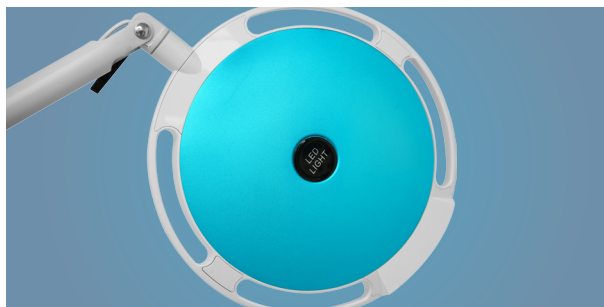

Light source components

The internal light source components are composed of ABS+aluminum substrate structure, offering excellent heat dissipation, lightweight to reduce failure rate.




Lamp head adjustment

Up and down 90°, horizontal flip 300°.



Height adjustment


Height is adjustable, can match the height of the patient easily.

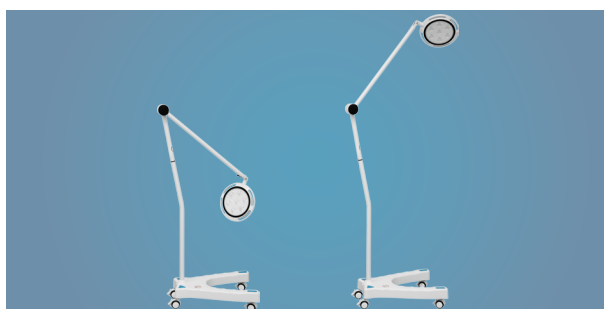



Base

U-sharp design, 4pcs 2-inch TPU universal silent wheels, easy and flexible to move.

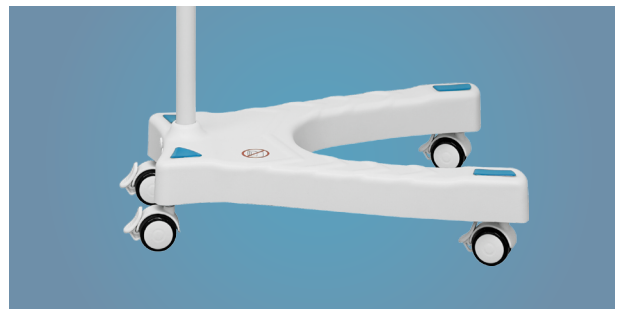

Control panel

Button operation is clear and straightforward.






CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited

Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC

17800 CASTLETON ST STE
180.
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

📞 +86 18021231901 · 🕒 24/7

✉️ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue