



SK-L306

WALL HANGING
MEDICAL LIGHT



SKSurgitech™
by Saikang Medical

Care you precisely™

SK-L306

FEATURES AND HIGHLIGHTS

Technical parameters

▪ Voltage	AC90-240V 50Hz~60Hz
▪ Bulb power	LED 3Wx12
▪ Color temperature	4000-5000K
▪ Bulb life	≥ 50000h
▪ Working distance	700mm
▪ Center illumination	≥45000Lux
▪ Minimum spot diameter	φ80±10mm
▪ Substrate temperature	≤55°C
▪ Lamp holder diameter	200mm
▪ Net weight	7.5KG
▪ Packing size	1040x335x230mm
▪ Packing material	carton

Description

This product adopts 12 pieces 3W LED bulbs as light source. High brightness, linear adjustable, daylight color, long lifetime. Base castor configuration, easy to move; Flexible tube bending, lamp holder 300° free rotation can get multi-angle lighting.



Application

General Surgery, Gynecological Examination, Ophthalmology, Dental, Cosmetology and Pet medical

Floodlighting

The maximum exposure area means the minimum shadow




Adjust the brightness

The brightness can be adjusted from maximum to off



Lamp holder axial

300° free rotation, multi-angle lighting

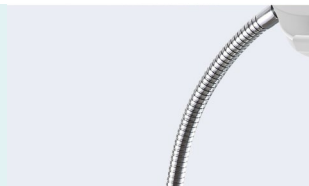

Castors in the back

Mobile convenient



LED light

Higher brightness, lower energy consumption and longer service time



Optional

Wall mounted type and Ceiling type





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi,
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai,
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎️ +86 18021231901 · 🕒 24/7
✉️ export@sksurgitech.com
🌐 www.sksurgitech.com



Company video



General catalogue