



Care you precisely!



SK-LLD04D

LED OPERATION LAMP

www.sksurgitech.com

SKSurgitech[™]
by Saikang Medica



SK-LLD04D

FEATURES AND HIGHLIGHTS

The LED surgical shadowless lamp is composed of multiple lamp bulbs. The bulb has a longer life than the traditional shadowless lamp. It is round or petal-shaped. It is fixed on the balance arm suspension system. It is stable in positioning and can be moved vertically or cyclically. It can meet different heights during surgery. The petal-shaped surgical shadowless lamp is also a kind of LED. It is introduced separately because the lamp cap is petal-shaped and has a special shape.



Technical parameters

▪ Diameter of lighthead (mm):	700
▪ Illumination at 1m (Lux):	160000
▪ Brightness control range:	10-100% stepless adjustment
▪ Color temperature (Kelvin):	3500-5000
▪ Color rendering index (CRI):	≥96
▪ Temperature rise (Surgeon head)(°C):	≤2
▪ Temperature rise (Surgical field)(°C):	≤1
▪ Total irradiance Ee (W/m2):	≤380
▪ Radiant power Ee/Ec (mW/m2lux):	≤4
▪ Effective focusing depth (cm):	100
▪ Lighting spot size (cm):	22-31
▪ L1+L2 depth of illumination (cm):	100
▪ Bulb type:	LED

▪ Bulb brand:	OSRAM
▪ Bulb quantity:	112
▪ Bulb power (w):	1
▪ Average bulb life (hours):	≥50000

Technical configuration

<input checked="" type="checkbox"/> Rotating arm	1 set
<input checked="" type="checkbox"/> Balance arm	1 set
<input checked="" type="checkbox"/> C-arm	1 set
<input checked="" type="checkbox"/> Handle	1 pcs

Material technology

The coating is sprayed with world-leading Akzo Nobel powder coating and anti-glare treatment, which meets the demanding requirements of surgical lighting. The design of the lamp holder focuses on reducing the turbulence. The streamlined lampshade of the whole lamp arm is carefully designed according to the principle of aerodynamics, which is fully in line with The standard of a modern new operating room.




Bulb

The bulbs are imported from Germany OSRAM brand, using a special optical system, 360 ° uniformly illuminated on the observed object, no shadows, true reproduction of the organization color, and achieve constant color temperature under any light intensity conditions. The material is hard and no mercury pollution. Also, the emitted light does not contain infrared and ultraviolet radiation pollution.




Color temperature&illumination

Using a new type of LED cold light source, the illuminance can reach 160,000 Lux, which can realize 16-level brightness adjustment; the color temperature is in the range of 3500-5000K, and four-level color temperature adjustment can be achieved. While meeting the high illuminance, adjust the lighting parameters according to the needs of different surgeons, so that the surgeon's perception of light is soft and no glint.

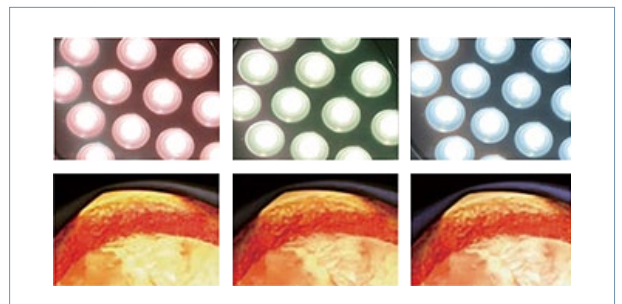


Color rendering index

The color rendering index reaches 96, which truly reflects the color of human tissues. It is suitable for various surgical scenes and greatly reduces the visual fatigue of medical staff due to long-term operations.

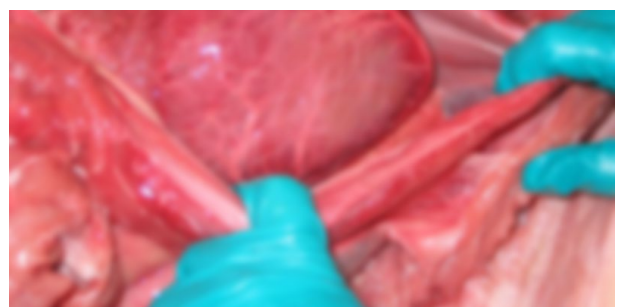

Temperature rise (surgical field)

LED cold light source has no infrared radiation, coupled with our nano-coated radiator, created our excellent heat dissipation effect. It solves the problem of the doctor's head temperature being too high during the operation, and creates a more comfortable operation environment for the doctor; it also effectively avoids the loss of water and infection of the surgical wound, which is conducive to the rapid recovery of the patient's wound after the operation.

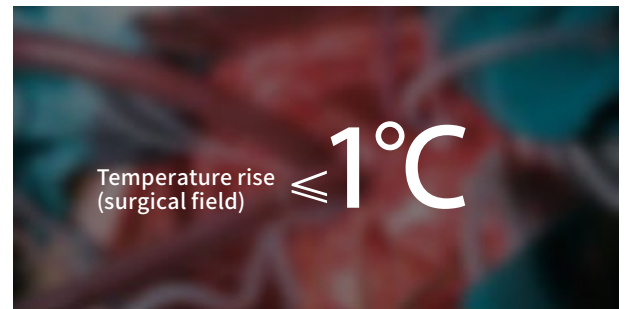


Electric focus

Compared with the traditional manual focusing function, this shadowless lamp is equipped with an electric focusing function, which is more convenient in operation.



Overhaul unit

Special overhaul unit, photoelectric separation design, is for quick repair, which keeps the circuit board away from the heat source. What's more, it increases the service life, and improve the convenience at the same time. The control box can also be easily disassembled for maintenance. The modular lamp body is made of high-quality aluminum alloy precision casting, and the light source can be replaced independently.

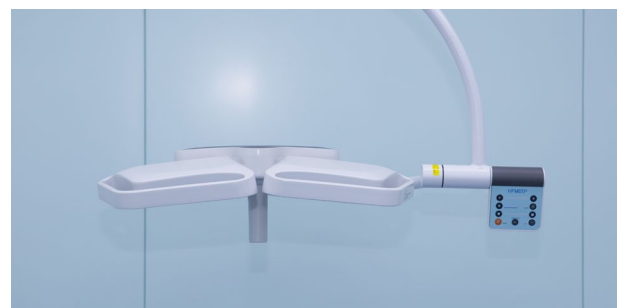

Operating handle


The detachable disinfection handle is made of high-strength ABS, which can be disinfected and sterilized at a high temperature above 135°, ensuring the durability and reliability of the handle.






LIGHT PATTERN

Intelligence of the adjustable lighting color temperature

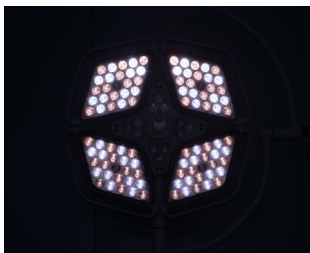


Normal mode

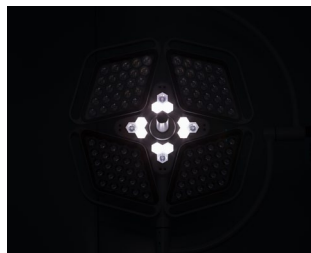


One-key cavity mirror mode

With one-key cavity mirror mode, red light compensation and green light compensation functions, and optional shadow compensation function.

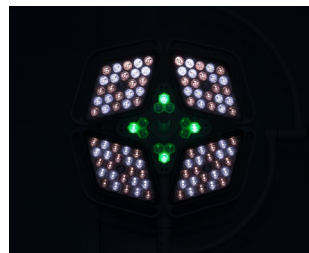

Green light mode

Suitable for surgery requiring high color rendering index, used for cardiothoracic surgery, neurosurgery, etc. Green light reduces the flare to the eyes caused by strong light.

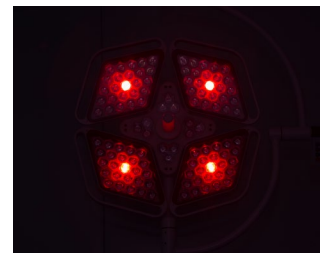

Red light mode

As a compensation function.



 **CONSTANT
INNOVATION**



Design



Research



Manufacturing



Maintenance



Service

 **WORLDWIDE
BRANCHES**

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa





Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi,
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech
by Saikang Medich

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 +86 18021231901 ·  24/7
 export@sksurgitech.com
 www.sksurgitech.com



Company video



General catalogue