


# SK-PT003

DOUBLE ARMS MEDICAL BRIDGE  
CEILING PENDANT

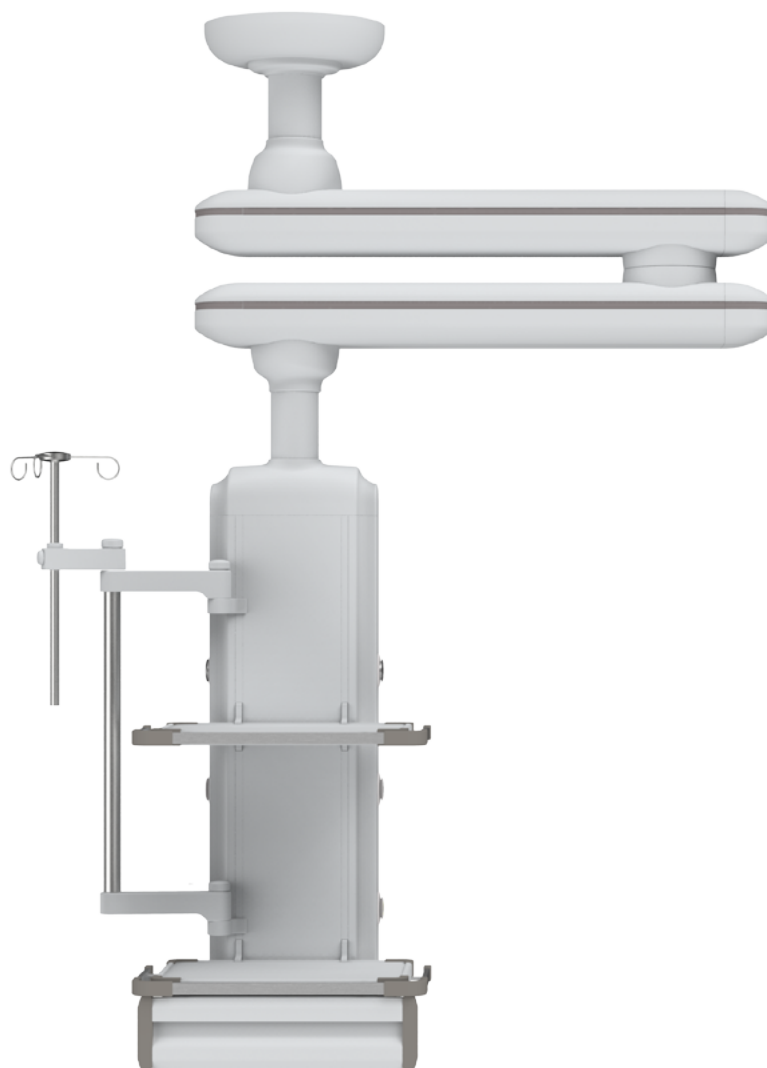

# SK-PT003

## FEATURES AND HIGHLIGHTS

The double arms medical bridge ceiling pendant is more flexible in moving and positioning, compared with the single one, longer treatment coverage radius, meeting the daily load-bearing requirements and using strong material, made of imported antibacterial plastic powder, ensuring the cleanliness of the operating environment. Its rich accessories matching various systems, greatly extending the clinical practicability.

Product Application: Suitable for OR room.

### Optional Accessories:

- \* Net basket - Stainless steel net basket
- \* We manufacturer and provide different standard of medical gas outlet: British Standard, American standard, German standard, French standard, Japanese Standard, etc.
- \* Different standard of medical gas outlets & adapters: British standard, French standard, German standard, Diss, Ohmeda, Chemetron, and etc.
- \* Inhalers - All kinds of oxygen inhalers (disposable, simple, external hanging, etc.)
- \* Negative pressure regulating valve
- \* Waste liquid collection device (external hanging type)

### Material and process

#### Appearance Design

The design of double hanging arms style, saves ground space; the whole bridge pendant is fully sealed, no sharp angle to prevent the impact.



#### Painting

Imported plastic powder, electrostatic painting process, painting thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



#### Welding

Panasonic Robot ensure 360° full smooth welding.



#### Material

The tensile strength of aluminum alloy frame is more than 295MPa, some connecting parts are made of cast iron  $\geq 200$ MPa, and the net load is more than 300kg; one-off extrusion molding, and the surface undergoes the primary oxidation treatment.



## Function

### Braking System

The pneumatic brake system, using compressed air to brake hanging arms and medical pendants, one press to release the brake.



### Suspension System

The rotation angle of medical pendant and hanging arm is 340°, built-in rotation limit device to prevent drift.



## Configuration

### Rotating Arm

The double rotating arms cover 1400mm treatment radius, can be used flexibly during treatment, meeting the surgical needs of all departments.



### Utility Model Patent

Negative pressure terminal obtained the Utility Model Patent CN202478404U.

专利内容由知识产权出版社提供

专利名称: 负压吸引气路终端  
 专利类型: 实用新型专利  
 发明人: 贾天升  
 申请号: CN202420028673.8  
 申请日: 2024/05/05  
 公开号: CN202478404U  
 公开日: 2024/10/10



摘要: 本实用新型公开一种负压吸引气路终端,包括两种结构,包括所述上盖设有负压吸引终端,在所述上盖的上部设有负压吸引终端,在所述上盖的上部设有负压吸引终端,在所述上盖的上部设有负压吸引终端。

### IV Pole with Four Hooks

The length of IV pole is 600mm, can be retractable, with four hooks, single hook can bear 4 kg.



### Drawer

Silent self-priming track, no damage after 10,000 times push-pull tests, and the handle position is integrally formed by pressing with mould.



### Instrument Platform

The size of the table top is 450\*500mm, and the height of table surface can be adjusted. The table top is made of aluminum alloy panel with iron square tube inside, can bear ≥50KG; with external aluminum alloy hanging frame.



### Medical Gas

It is equipped with medical gas on both sides of the cabinet and column. The following configurations can be customized,  
 1. The number of medical gas  
 2. Interface standard and position  
 Gas port adopts secondary seal, three states (on, off, pull), support maintenance operations while it is working.



### Equipotential Socket

Reduce the potential difference between indirect contact voltage and different metal parts.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@sksurgitech.com

🌐 www.sksurgitech.com



Company video



General catalogue