


SK-PT005

ANESTHESIA MEDICAL BRIDGE
CEILING PENDANT

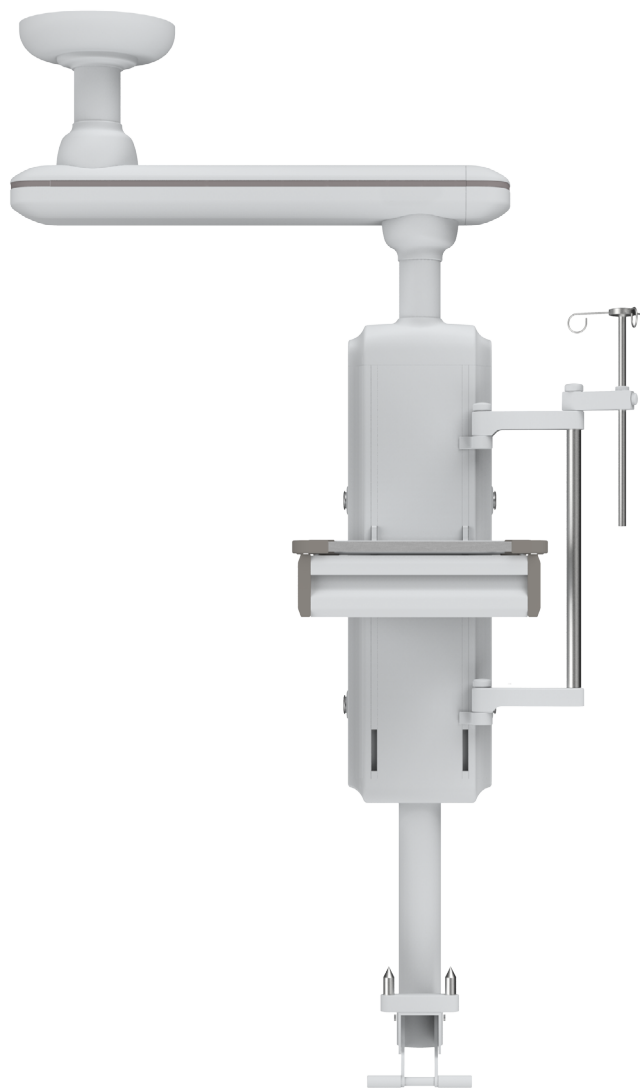

SK-PT005

FEATURES AND HIGHLIGHTS

The anesthesia medical bridge ceiling pendant is equipped with hanging holder for anesthesia machine, bearing 250 kg of anesthesia machine, with international AGSS jet exhaust gas, suitable for anesthesia surgery.

Product Application: Suitable for OR room.


Optional Accessories:

- * Net basket - Stainless steel net basket
- * We manufacturer and provide different standard of medical gas outlet: British Standard, American standard, German standard, French standard, Japanese Standard, etc.
- * Different standard of medical gas outlets & adapters: British standard, French standard, German standard, Diss, Ohmeda, Chemetron, and etc.
- * Inhalers - All kinds of oxygen inhalers (disposable, simple, external hanging, etc.)
- * Negative pressure regulating valve
- * Waste liquid collection device (external hanging type)

Material and process

Appearance Design

The design of single hanging arm style, saves ground space; the whole bridge pendant is fully sealed, no sharp angle to prevent the impact.

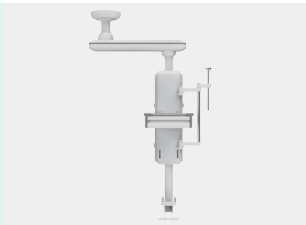

Painting

Imported plastic powder, electrostatic painting process, painting thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



Welding

Panasonic Robot ensure 360° full smooth welding.



Material

The tensile strength of aluminum alloy frame is more than 295MPa, some connecting parts are made of cast iron ≥ 200 MPa, and the net load is more than 300kg; one-off extrusion molding, and the surface undergoes the primary oxidation treatment.



Function

Braking System

The pneumatic brake system, using compressed air to brake hanging arms and medical pendants, one press to release the brake.



Lifting System

The function column built-in Jiechang motor of model JC35D1, with SAA/TUV/PSE/CE/UL/KC certificates, ensure safety. Motor noise under 50dB, lifetime 20,000 times, electric lifting distance 300mm; accurate positioning without drift, meeting requirements of different nursing.



Suspension System

Medical pendant and hanging arm, rotation angle 340°, built-in rotation limit device to prevent drift.



Configuration

IV Pole with Four Hooks

The length of IV pole is 600mm, can be retractable, with four hooks, single hook can bear 4 kg.



Drawer

Silent self-priming track, no damage after 10,000 times push-pull tests, and the handle position is integrally formed by pressing with mould.



Instrument Platform

The size of the table top is 450*500mm, and the height of table surface can be adjusted. The table top is made of aluminum alloy panel with iron square tube inside, can bear ≥ 50 KG; with external aluminum alloy hanging frame.



Medical Gas

It is equipped with medical gas on both sides of the cabinet and column. The following configurations can be customized,
1. The number of medical gas
2. Interface standard and position
Gas port adopts secondary seal, three states (on, off, pull), support maintenance operations while it is working.



Rotating Arm

The single rotating arm can cover 1000mm treatment radius, can be used flexibly during treatment, meeting the surgical needs of all departments.



Anesthesia Machine Holder

Anesthesia machine holder, bearing 250 kg of anesthesia machine.



Special Gas Terminal

With carbon dioxide gas terminal for anesthesia surgery, laughing gas terminal and Jet-discard AGSS collection terminals.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai,
United Arab Emirates.
Tel.: +20 128 8962497

Healthward^{uk}

Chengdu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No.1, 5th Floor, Building 5, No.360, Wuke West 5th Road ·
Wuhou District (Chengdu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@health-ward.com

🌐 www.health-ward.com



Company video



General catalogue