


# SK-PT006

## SUPPLY COLUMN

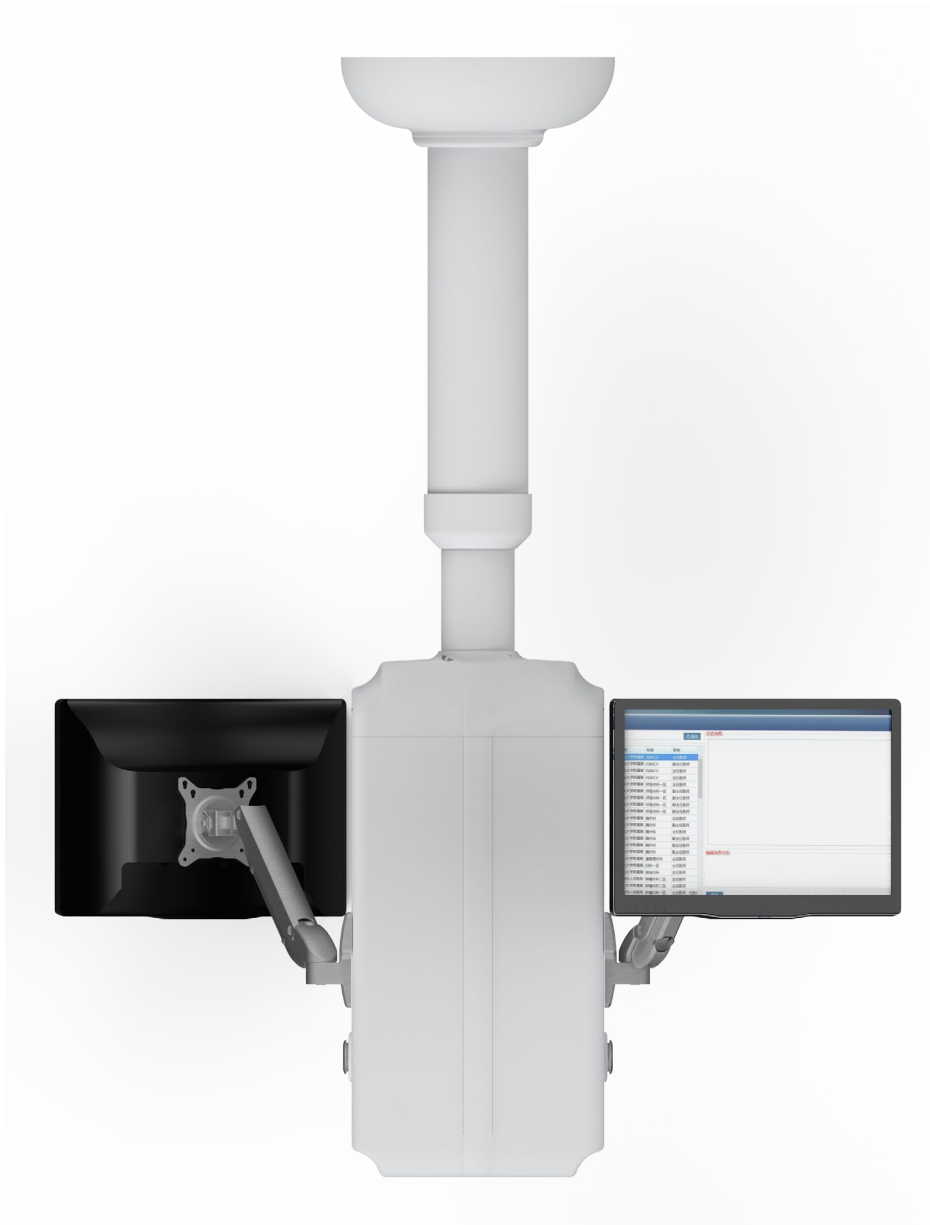

# SK-PT006

## FEATURES AND HIGHLIGHTS

Working table supports can be placed on both sides of the supply column, one-stop operation for medical staff, suitable for supply rooms.

Product Application: Suitable for supply room.

### Optional Accessories:

- \* Net basket - Stainless steel net basket
- \* We manufacture and provide different standard of medical gas outlet: British Standard, American standard, German standard, French standard, Japanese Standard, etc.
- \* Different standard of medical gas outlets & adapters: British standard, French standard, German standard, Diss, Ohmeda, Chemetron, and etc.
- \* Inhalers - All kinds of oxygen inhalers (disposable, simple, external hanging, etc.)
- \* Negative pressure regulating valve
- \* Waste liquid collection device (external hanging type)

#### Appearance Design

The design of ceiling type function column, saves ground space; the whole function column is fully sealed, no sharp angle to prevent the impact.



#### Painting

Imported plastic powder, electrostatic painting process, painting thickness 0.12mm, brightness 60°, paint can resist 50kg impact.




#### Welding

Panasonic Robot ensure 360° full smooth welding.



#### Material

The tensile strength of aluminum alloy frame is more than 295MPa, some connecting parts are made of cast iron  $\geq 200$ MPa, and the net load is more than 300kg; one-off extrusion molding, and the surface undergoes the primary oxidation treatment.

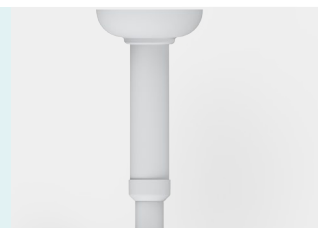

#### Suspension System

Medical pendant rotation angle 340°, built-in rotation limit device to prevent drift.




#### Working Table Support

Equipped with 2 sets of 349° rotating 300mm diameter working table supports.




#### Equipotential Socket


Reduce the potential difference between indirect contact voltage and different metal parts.

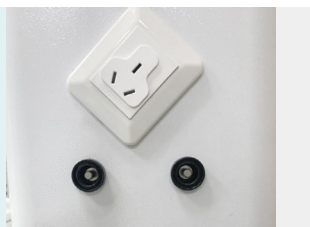

#### Medical Gas

It is equipped with medical gas on both sides of the cabinet and column. The following configurations can be customized,  
1. The number of medical gas  
2. Interface standard and position  
Gas port adopts secondary seal, three states (on, off, pull), support maintenance operations while it is working.

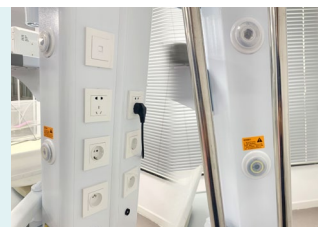

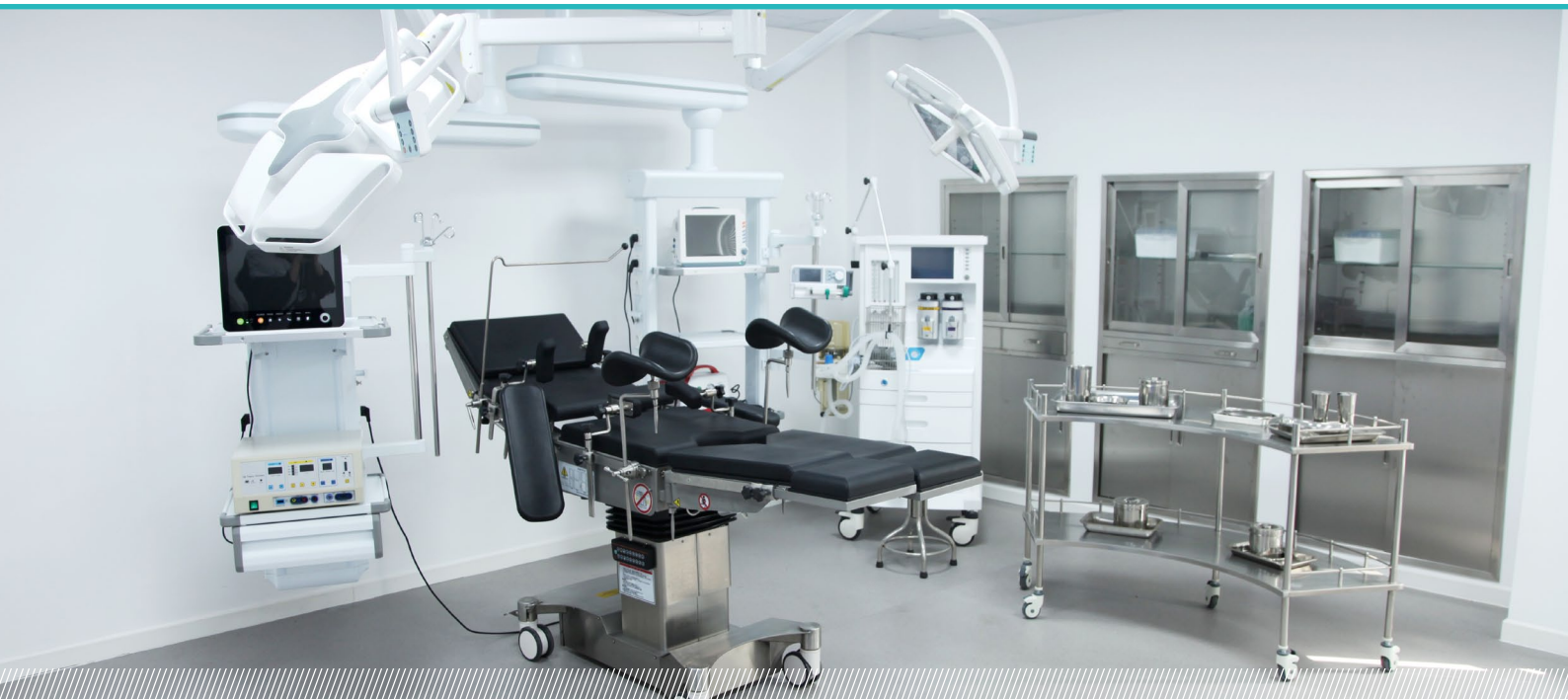

## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi,  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai,  
United Arab Emirates.  
Tel.: +20 128 8962497

# Healthward

Chengdu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No.1, 5th Floor, Building 5, No.360, Wuke West 5th Road ·  
Wuhou District (Chengdu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@health-ward.com

🌐 www.health-ward.com



Company video



General catalogue