


# SK-SY05

## OSCILLATOR

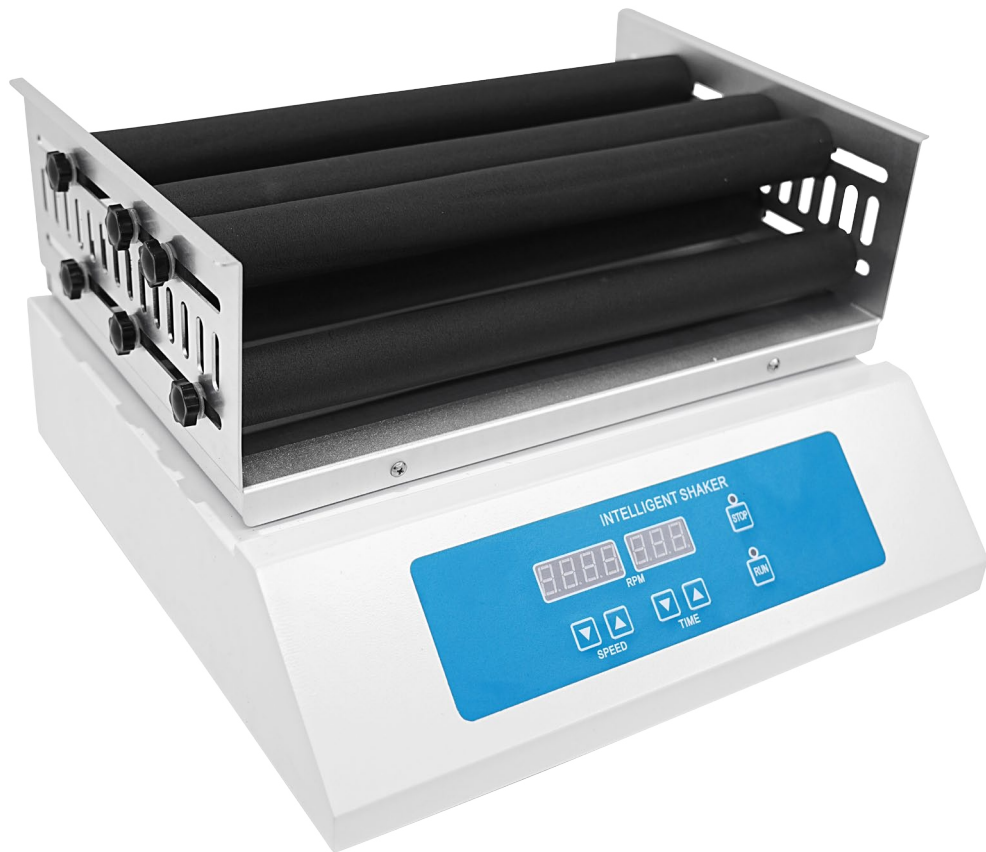

# SK-SY05


## FEATURES AND HIGHLIGHTS

**Design:** The product is suitable for laboratories in scientific research departments such as medical and health, pharmaceutical industry, biochemistry, industrial production, and agricultural science. It is an experimental instrument used for swirling and mixing various samples under the same environmental conditions, also known as a rotary oscillation mixer.

**Working principle:** A permanent magnet DC motor drives a rotating tray to oscillate and mix experimental samples.

### Technical parameters

▪ Oscillation mode	rotatory-type
▪ Work table size (mm)	420x290
▪ Speed adjustable range	start-300 rpm
▪ Timing	0~990 minutes
▪ Amplitude	20mm
▪ Power supply	110v 60hz/220v 50hz
▪ Overall dimensions (mm)	420x410x260
▪ Gross/net weight	30/18kg

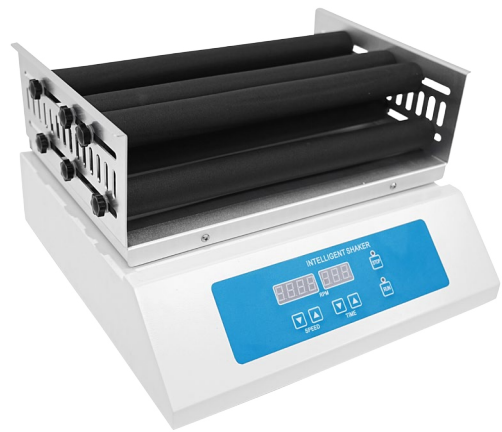


#### Operation panel

LCD display, button operation is simple.



#### Material

Cold-rolled steel material, electrostatic spray coating, no paint peeling, corrosion resistant.

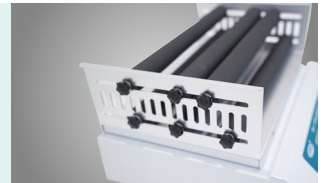

#### Clamp-style tray

Clamp-style tray, adjustable distance.



#### Motor

The permanent magnet DC motor is driven by a precision electronic speed control circuit, and it has the characteristics of smooth rotation, low noise, and adjustable speed.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi,  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@sksurgitech.com

🌐 www.sksurgitech.com



Company video



General catalogue