



SK-SY16

MAGNETIC HEATING STIRRER



SKSurgitech™
by Saikang Medical

Care you precisely

SK-SY16


FEATURES AND HIGHLIGHTS

Design: Instruments are widely used in industries such as medicine and health, clinical testing, biochemistry, and environmental protection. They are one of the commonly used analytical instruments in physical and chemical laboratories. They are used for liquid mixing, mainly for stirring or heating and stirring low-viscosity liquids or solid-liquid mixtures.

Working principle: By using the principle of magnetic field's like poles repel and opposite poles attract, a magnetic stirrer placed in a container is driven in a circular motion through the magnetic field to stir the liquid. Together with a temperature control system, the liquid can be heated to the desired experimental temperature.

Technical parameters

▪ Speed range	start ~1600 rpm
▪ Stirring power	25w
▪ Temperature range	rt+5 ~99°C
▪ Heating power	120w
▪ Table size(mm)	φ120
▪ Stirring capacity	100~1000ml
▪ Power supply	110v 60hz/220v 50hz
▪ Machine size(mm)	225x160x130
▪ Gross/net weightt	2.5/2Kg

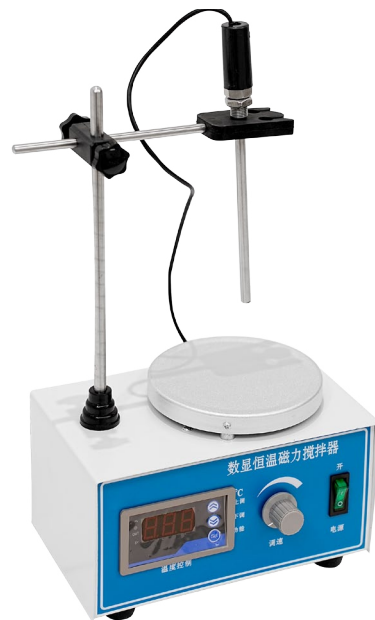

Sample tray

The stainless steel sample tray is sturdy, corrosion-resistant, and easy to clean.




Display screen

The LCD liquid crystal display screen shows clear and beautiful content, and the reading is accurate and intuitive.



Shell material

The box body is made of thick cold-rolled steel plate, which is firm and reliable. The electrostatic spray-coated shell does not peel off paint and is corrosion-resistant.

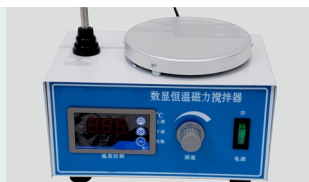

Pole

Stainless steel uprights are used, which are sturdy, durable, and corrosion-resistant.



Temperature Sensor

An external temperature sensor is used to monitor and control the sample temperature in real-time.



Strong magnet

Strong magnetic magnets are used, which have a significant stirring effect.





CONSTANT INNOVATION




Design



Research



Manufacturing

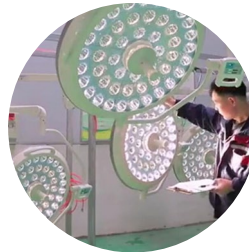

Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi,
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai,
United Arab Emirates.
Tel.: +20 128 8962497

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ export@sksurgitech.com

🌐 www.sksurgitech.com



Company video



General catalogue