



SKR-R01

INFORMATION NURSING CART

Registration code: SKM-A



SKSurgitech™

by Saikang Medical



SKR-R01

FEATURES AND HIGHLIGHTS

Technical Parameters

▪ Maximum external dimension W×D×H	520×620×1315mm
▪ Table size W×D	480×60mm
▪ Table height from the ground	935mm
▪ Weight of the whole machine	50KG (without computer)
▪ External dimensions of the box W×D×H	480×460×535mm
▪ 1st drawer inner dimension W×D×H	390×380×40mm
▪ 2nd, 3rd and 4th drawer inner dimension W×D×H	390×380×60mm
▪ Pull-out sub-table size L×W	260×260mm
▪ Display mounting hole	75mm×75mm, 100mm×100mm
▪ Constant force lifting stroke	90mm

Optional Accessories

- Hand sanitizer basket
- Mouse hanger
- ABS storage basket
- Trash can
- Sharps box
- Power plug hanger
- Probe cup



Material

The table top is made of ABS medical grade plastic, the four-layer all-aluminum box is made of high-voltage electrostatic spraying process, and imported environmentally friendly antibacterial plastic powder is used to ensure that the product meets the medical and health requirements, and the surface is matte and non-dazzling. And equipped with an integrated die-cast aluminum handle.



Upper Column

The detachable upper column can carry a monitor or an all-in-one computer up to 10kg, and the universal support function can complete 180° left and right rotation and 45° up and down inclination. The up and down stroke of the column is as high as 9cm, with a constant force lifting adjustment for easy operation, and various equipment threading can be hidden in the column to make the appearance more refreshing.



Table Top

The ABS medical-grade plastic countertop is easy to clean and anti-corrosion, without any spraying on the surface, which will not affect the medicine. The grooves are set around to prevent the dripping of medical liquid, and there is a countertop enclosure to prevent items from falling.




Sub Table Top

The pull-out sub-table provides the area expansion of the desktop, which can be temporarily convenient for work.



Keyboard Drawer

Hidden keyboard drawer that can be stored inside the countertop when not in use.

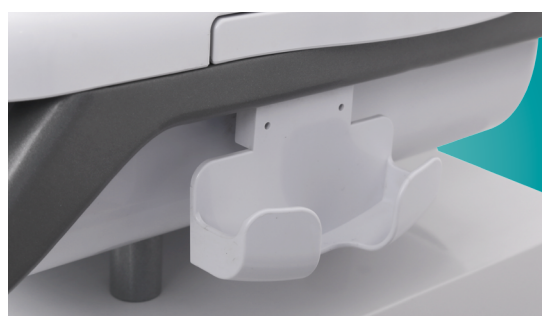

Electromechanical Layer

The distance between the electromechanical layer and the medicine box is greater than 5cm, which facilitates the heat dissipation of the electromechanical equipment and protects the safety of medicines.



Drawer

Four-layer large space self-priming drawer, the inner layer of the drawer is equipped with a plastic insert board, which can be combined in any size to store medicines of different specifications.



Handle

The ring-shaped anti-collision handle, surrounded by solid metal handles with rounded corners, is located at the outermost periphery of the whole vehicle, which can not only injure other personnel in the event of a collision, but also protect the safety of the stroller.



Caster

4-inch dual-wheel medical grade silent wheels, all with brakes, after tens of thousands of kilometers of long-term durability tests, flexible and durable.





CONSTANT INNOVATION

Design

Research

Manufacturing

Maintenance

Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd

Sala Thammasop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited

DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya Naberezhnaya,
Russia
TEL: +7 (919)300-00-06

Mid East Branch


Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SKSurgitech™
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

 +86 18021231901 ·  24/7

 export@sksurgitech.com

 www.sksurgitech.com



Company video



General catalogue