



# SKR-R02

## MOBILE MULTI-FUNCTION INFORMATION SEARCH CART

Registration code: SKM-A



**SKSurgitech™**

by Saikang Medical

PRODUCT  
DATA SHEET

# SKR-R02

## FEATURES AND HIGHLIGHTS

### Technical Parameters

▪ Dimensions L×W×H	620×540×1460mm
▪ Work table L×W	460×480mm
▪ Weight of the whole machine	45KG (without computer)
▪ Sub-table size L×W	260×260mm
▪ Trolley lift stroke	200mm
▪ Display lifting stroke	90mm
▪ Caster design	4 inch medical grade silent wheel
▪ Cabinet parameters	two layers
▪ Drawer dimensions	480×340×295mm
▪ Layer drawer inner size	390×275×55mm
▪ The inner size of the drawer on the second floor	390×275×140mm

### Optional Accessories

- Hand sanitizer box
- Mouse hanger
- Sundries box
- Trash can



## Material

The table top is made of ABS medical grade plastic, the two-layer all-aluminum box is made of high-voltage electrostatic spraying process, and imported environmentally friendly antibacterial plastic powder is used to ensure that the product meets the medical and health requirements, and the surface is matte and non-dazzling. And equipped with an integrated die-casting aluminum handle and base, The column material is made of medical grade aviation aluminum alloy.



## Upper Column

The detachable upper column can carry a monitor or an all-in-one computer up to 10kg, and the universal support function can complete 180° left and right rotation and 45° up and down inclination. The up and down travel of the column is as high as 15cm, with a constant force lifting adjustment for easy operation, and various equipment threading can be hidden in the column to make the appearance more refreshing.



## Table Top

The ABS medical-grade plastic countertop is easy to clean and anti-corrosion, without any spraying on the surface, which will not affect the medicine. The grooves are set around to prevent the dripping of medical liquid, and there is a countertop enclosure to prevent items from falling.




## Sub Table Top

The pull-out sub-table provides the area expansion of the desktop, which can be temporarily convenient for work.



## Keyboard Drawer

Hidden keyboard drawer that can be stored inside the countertop when not in use.

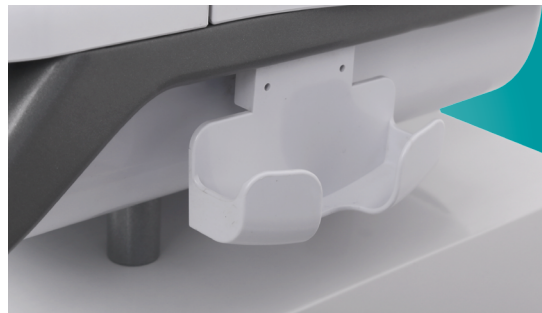

## Electromechanical Layer

The distance between the electromechanical layer and the medicine box is greater than 5cm, which facilitates the heat dissipation of the electromechanical equipment and protects the safety of medicines.



## Drawer

Two-layer large space self-priming drawer, the inner layer of the drawer is equipped with a plastic insert board, which can be combined arbitrarily to store medicines of different specifications.



## Handle

The ring-shaped anti-collision handle, surrounded by solid metal handles with rounded corners, is located at the outermost periphery of the whole vehicle, which can not only injure other personnel in the event of a collision, but also protect the safety of the stroller.



## Lower Column

The column is equipped with a professionally adjusted air rod, which is easy to rise and fall, and can be quickly adjusted to a suitable height, adapting to different sitting and standing postures, with a stroke of up to 20cm. And optional electric lifting system.



## Chassis

All-metal stretched X-shaped chassis, one-time cast aluminum molding, to ensure the center of gravity and level of the vehicle.



## Casters

4-inch dual-wheel medical grade silent wheels, all with brakes, after tens of thousands of kilometers of long-term durability tests, flexible and durable.



## Monitor (optional)

Split computer: 21.5 inches or 23.8 inches  
All-in-one computer: 23.8 inches and above





## CONSTANT INNOVATION

Design

Research

Manufacturing

Maintenance

Service



## WORLDWIDE BRANCHES

### South America Branch

#### Saikang Medical Technology, Inc

W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

#### Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

#### Saikang Medical Technology Ltd

Sala Thammasop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

#### Uk Healthward International Limited

DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 2021111106

### Russia Branch

#### LLC 'SAIKANG'

454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

#### Saikang Medical Supplies (Mid East) LLC

Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

**SKSurgitech™**  
by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

 +86 18021231901 ·  24/7

 export@sksurgitech.com

 www.sksurgitech.com



Company video



General catalogue