



SKR-R04

WARD ROUNDS CART

Registration code: SKM-A



SKSurgitech™
by Saikang Medical



PRODUCT
DATA SHEET

SKR-R04

FEATURES AND HIGHLIGHTS

Product Introduction

The mobile doctor workstation based on wireless network integrates highly integrated computer, wireless communication, intelligent power management, mobile management technology, and is suitable for a variety of clinical work environments; whether it is outpatient, emergency, inpatient ward, operating room or ICU, it can help clinical Medical workers realize intelligent mobile ward rounds.

Technical Parameters

▪ Dimensions L×W×H	530×520×730mm (table from the ground)
▪ Work table L×W	420×380mm
▪ Overall weight	17.2kg
▪ Lifting method	table lifting size 30mm
▪ Caster design	pneumatic rod
▪ Lifting size	table lifting size 200mm
▪ Caster design	4 inch medical grade silent wheel
▪ Monitor arm without constant force lifting height	410mm
▪ Monitor arm without constant force lifting bearing	20kg



Material

The table top and column are made of aviation-grade aluminum alloy, the base is made of one-time molding all-cast aluminum material, the whole car adopts high-voltage electrostatic spraying process, and imported environmentally friendly antibacterial plastic powder is used to ensure that the product meets the medical and health requirements, and the surface is matte and non-dazzling.



Upper Column

The detachable upper column can carry an all-in-one computer up to 10kg, and the universal support function can complete 180° left and right rotation and 45° up and down inclination. The up and down travel of the column is as high as 15cm, with a constant force lifting adjustment for easy operation, and various equipment threading can be hidden in the column to make the appearance more refreshing.



Table Top

Aviation grade aluminum alloy table, easy to clean and anti-corrosion, with partition design on both sides to prevent items from falling.



Lower Tray

The lower tray can be equipped with a printer, which can simply print the patient's medical records and some materials while moving the information-based ward rounds.



Lower Column

The column is equipped with a professionally adjusted air rod, which is easy to rise and fall, and can be quickly adjusted to a suitable height, adapting to different sitting and standing postures, with a stroke of up to 20cm.




Chassis

All-metal stretched X-shaped chassis, one-time cast aluminum molding, to ensure the center of gravity and level of the vehicle.



Casters

4-inch dual-wheel medical grade silent wheels, all with brakes, after tens of thousands of kilometers of long-term durability tests, flexible and durable.

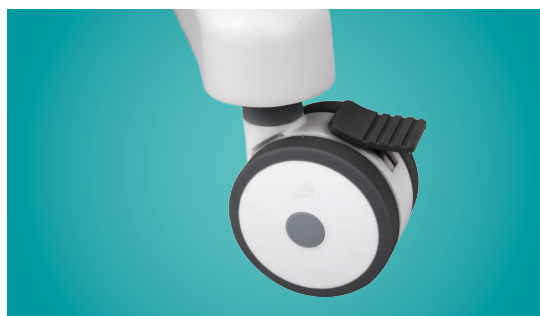



CONSTANT INNOVATION

Design

Research

Manufacturing

Maintenance

Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc

W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd

Sala Thammasop Subdistrict
Thawi
Watthana District
Bangkok Thailand
TEL: +66 642741046

UK Branch

Uk Healthward International Limited

DEPT 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +31 (0) 2021111106

Russia Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya Naberezhnaya,
Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC

Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SKSurgitech™

by Saikang Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

+86 18021231901 · 24/7

export@sksurgitech.com

www.sksurgitech.com



Company video



General catalogue