

## NEWS




## CE Certification

NEW CE regulation MDR EU 2017745 in 2020

- Year-end Strategy Conference 4
- Thanks for your trust! 5
- Innovating: Y8t bed 6



Jiangsu Saikang Medical Equipment Co., Ltd has completed the registration in the Netherlands in 2020 for six categories of Electric Hospital Bed, Manual Hospital Bed, Electric Operating Table, Operating Table, Medical Trolley and Stretcher, and successfully passed the certification of the new CE regulation MDR EU 2017/745, which ensures the free circulation of our products in the EU market. Our product registration information can be found on official website of EU authorized Representative: Shanghai Sungo Medical Technology Co., Ltd. through the company name and certification document number.

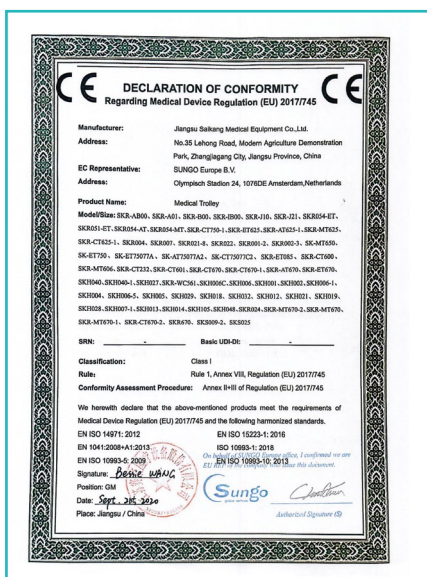

# Year-end Strategy Conference



On January 16, 2021, Saikang held the conference to state results we achieved in 2020 and made more strategic decisions to meet possibilities in the coming 2021. General manager Michael and senior leaders from Saikang different departments participated in this conference.

In the past year, electric beds, manual beds, and all our complementary products have been of great help to the 95 countries around the world. Normal operations of many large orders can't do without active cooperation of all departments. We

show sincere thanks to them for their great contributions.

Now, all departments will be well prepared for possible challenges in 2021 and look forward to bright future.



# Thanks for **your trust!**




## Shew La Min Hospital

Shwe La Min Hospital is a dedicated health care institution facilitating patients treatment with specialized medical and nursing staff and medical equipment. It was established in March 2009 to provide unsurpassed health and wellness services. The main objective is to achieve patients' satisfaction by providing premium quality health care. Shwe La Min Hospital equipped electric hospital beds from us in the wards mainly.



## LinnThitSar Hospital

We thank LinnThitSar Hospital for sharing these beautiful photos of Saikang products in this wonderful hospital of Myanmar. In order to improve health care service, we support complete product range including hospital beds with bed mattresses, transportation trolleys, overbed tables, night tables, accompany chairs, medical cabinets and anesthesia trolleys.

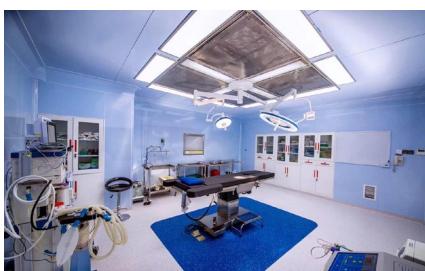

Oasis Series


# Y8t

## MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C


Size (LxWxH): 2140x1060x(480-770) mm

Product video

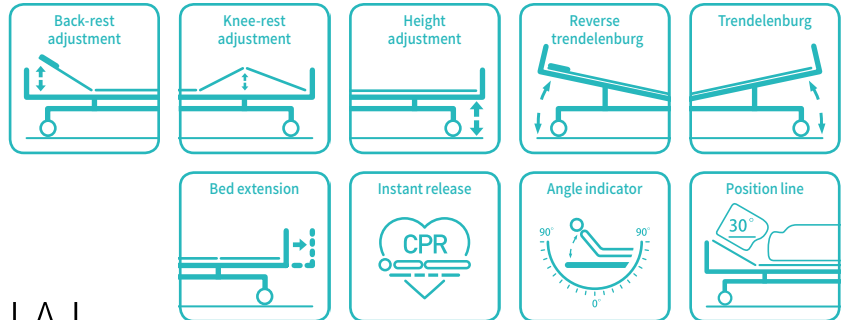

**Great multi-care**

Respects easy patient mobilisation: One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing.

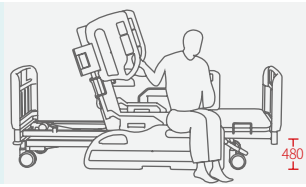

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding. Solid bed base made of 5.0 mm stamped rotating plate.



**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.



**Electric motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.



**Safety**

IEC60601 Safty Standard. One-button emergency stop switch can quickly stop all going movement.



**CPR**

CPR emergency reset handle to immediately flatten the bed back part; located on both sides of the bed. Electric CPR button is setted inside control panel, one button to make the back plate flat.



**Wheels**

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts. Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.



**Brake**

Three position braking system, transporting can maintain a straight line.



**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.



**Angle indicator**

Patent angle indicator, clearly show back-rest and tredelenbug position.




**Head & foot board**

Detachable PP head&foot board(Y19), with medical record card on foot panel, patent design, using the eccentric locking structure, one step quick removable.




**Siderail**

4 pcs PP siderail (T20), with patent angle indicators on each one, clearly show back-rest and tredelenburg position.

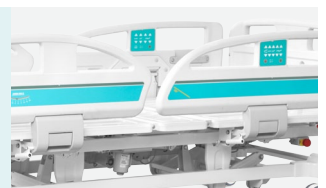


**Siderail controller**

Controls in the siderails: Numerous studies have shown that the built-in controller located in the head guardrail is easier to reach, allowing the patient to adjust the basic position of the bed through the controller. It is different between inside and outside, more safe and scientific.

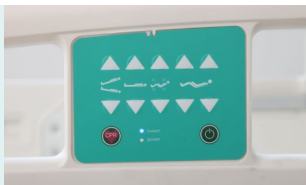

**Wall bumper wheel**

Φ70mm bumper wheels at the outermost periphery of the bed to protect the bed and wall, anti-collision and durable.

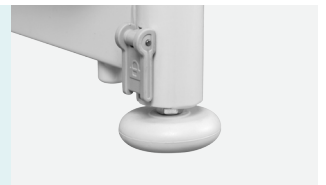

**IV pole**

Stainless Steel Column with plastic hook.



**Drainage hooks**

Mobile plastic drainage hooks.



# In brief

## Happy 2020 Christmas

On the afternoon of December 25th, all members of international sales and marketing department participated in Christmas celebration activity held by company. It was a good chance for all to enjoy a very wonderful movie together at the cinema. When night coming, along with members of several departments, we played games, enjoyed delicious food, sang and entered into the spirit of Christmas. It was a perfect Christmas.



## New Season Factory Live Show

We have started to do factory live show since last summer so as to be seen by more friends around the world. It takes time preparing promotion products and sparing work times. But it is worth being seen by the world through the ordinary camera. After all, we devoted a lot very much. Now, we have new live shows every month in 2021. Our professional team will show you around SAIKANG factory for a new look. Welcome to visit SAIKANG online.



## Stock Hospital Beds during Chinese New Year

Professional in providing medical solutions for medical institutions and health-care centers for nearly 20 years, now we have branches in USA, Nairobi, London, UAE, and Russia for good service. During Chinese new year holiday in 2021, we stock hundreds of hospital beds in advance in case of emergency health-care problems. Besides, our salesmen still serve customers anytime during holiday time. Product department have only 2 days off to keep 60% capacity to run. All we do is to fulfill customers' requirements.

**SAIKANG Medical**  
LEADING BRAND OF CHINA MEDICAL FURNITURE

**READY STOCK**  
Stock Code: 870098

**A1k**  
Manual Bed  
51K-CD1

75° & 110°  
Reclining



**Technical Parameters**

External Length	2110mm
External Width	960mm
Bed Height	500mm

**Technical Configuration**

<input type="checkbox"/> Manual Crank	1pc	<input type="checkbox"/> W/Pain Position	3pcs
<input type="checkbox"/> 3° Covered Castors	4pcs	<input type="checkbox"/> Drainage Hook	4pcs
<input type="checkbox"/> Adjustable Side Rail	1set	<input type="checkbox"/> ABS Anticollision	4pcs

**SAIKANG Medical**

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

📞 +86 18021231901 · 🕒 24/7

✉️ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video    General catalogue