

CELEBRATING
ANNIVERSARY

24th

· 2002-2026 ·



Healthier your life

Smart Mobile Medical Workstation Cart Solution

www.saikangmedical.com

SAIKANG Medical®

- Enhancing Workflow
- Improving Efficiency
- Empowering Healthcare Providers

STRENGTH OF SAIKANG

24

24 years focusing on design, production, service and technical consultation of medical care equipment

150

Located in Zhangjiagang City, 150km nearby Shanghai Economic Center

1000000

Saikang delivered 1,000,000 sets of products globally. Products all over Asia, Middle East, North America and Africa

400

There are more than 400 employees worldwide, more than 50 qualified engineers, and more than 100 trained technicians

245

Equiped with 245 sets of various automatic production facilities, Saikang follow strict standard control system on every production process, enable the traceable of each part

50000

The finished phase I production area covers more than 30,000 square meters, and phase II will be more than 20,000 square meters of advanced production area

31

Saikang offer comprehensive medical solutions with 31 product categories

562

Safe, efficient and easy-to-use medical solution including emergency, nursing and ICU wards of more than 500 series of products

32

With CE certification, ISO9001/13485 quality system certification, U.S. FDA and other 32 certifications

189

Ownership of 189 patents, continuously leading new standards in defining product quality and service efficiency

150

Since the establishment of the company, products and solutions have been used in more than 150 countries

1000

Saikang products are committed to creating a modern nursing environment, accumulatively serving more than 1,000 medical and health care organizations



LEADING BRAND OF MEDICAL FURNITURE

Patients in medical rehabilitation often have to face huge challenges. Many people who suffering from diseases or accidents are faced with retraining their minds and bodies and must adapt to the new environment. In this case, the positive atmosphere would offer decisive help. Hospital beds and furniture should provide the best support for patients and nurses. Considering this case, Saikang devoted to provide medical care solutions and perfect choices for medical institutions and rehabilitation centers, to help patients achieve rapid or long-term recovery.

Saikang was founded in 2002, now is a leading global provider of medical furniture and solutions. Head office locates in Suzhou, and branches in US, Dubai, Russia, Brazil, Thailand and Kenya, having 400 employees worldwide. Saikang has formed a large global network of research and development, marketing and service. Currently, Saikang products and medical solutions have been used in more than 149 countries and regions in the world.

Frontline Challenges in Daily Clinical Work



- Nurses and clinicians spend a significant amount of time walking back and forth between patient beds and nursing stations, reducing time available for direct patient care and lowering overall workflow efficiency.
- Paper-based documentation and fragmented digital systems delay real-time data entry, increase the risk of transcription errors, and make it difficult to access up-to-date patient information at the point of care.
- Traditional computer carts are often bulky and poorly balanced, making them difficult to maneuver in crowded clinical environments. Limited height adjustment and non-ergonomic design may lead to fatigue and long-term musculoskeletal strain.
- Strict infection-control standards require medical equipment with smooth, seamless, and easy-to-clean surfaces. Conventional carts with exposed cables, joints, and gaps increase cleaning time and cross-contamination risk.

These challenges highlight the urgent need for smarter, mobile, and ergonomically designed medical workstation carts that support efficient, safe, and compliant clinical workflows.

Why Choose Our Medical Workstation Carts?

True Mobility & Flexibility

The ultra-lightweight structure combined with medical-grade mobility components allows caregivers to move smoothly and effortlessly across wards, corridors, and patient rooms.

Equipped with 4-inch double-wheel medical-grade silent casters, including two lockable brake casters, the cart ensures both excellent maneuverability and reliable stability when stationary.

This design minimizes physical effort, reduces noise in clinical environments, and enables caregivers to focus more on patient care rather than equipment handling.



Ergonomic Design for Comfort and Safety

The column is equipped with a professionally calibrated gas spring, allowing smooth, effortless lifting and lowering with minimal force. Caregivers can quickly adjust the workstation to an optimal height for both sitting and standing positions, improving comfort during long shifts.

Additionally, the circular bumper handle made of sturdy metal with rounded edges is positioned around the outer perimeter of the cart, which helps prevent injury to surrounding personnel in the event of collisions while also protecting the cart itself from accidental damage.



Modular & Highly Customizable

The modular design allows the cart to be easily configured according to specific clinical workflows and departmental needs.

It supports a wide range of accessories, including printers, barcode scanners, drawers, medication boxes, CPU holders, and storage baskets, enabling seamless integration into nursing, medication administration, and documentation tasks.

This flexibility ensures that one platform can adapt to multiple clinical scenarios without the need for additional equipment.



Infection-Control Design

Designed with infection control as a priority, the cart features smooth, seamless surfaces that minimize dirt accumulation and simplify cleaning.

The antibacterial coating and cable-free layout reduce potential contamination points, making routine disinfection faster and more effective.

This design supports hospitals in meeting strict infection-control and hygiene standards, particularly in high-risk environments such as ICUs and operating areas.



TECHNICAL SPECIFICATIONS



Model	SKR-R01	SKR-R02	SKR-R04	SKR-R07
Length	500±10mm	530±10mm	535mm (±10mm)	530mm (±10mm)
Width	620±10mm	620±10mm	520mm (±10mm)	620mm (±10mm)
Height	1015±10mm	945±10mm	725-1045mm (±10mm)	855-1100mm (±10mm)
Tabletop Size	480*460±10mm	480*460±10mm	450×505mm (±10mm)	480×460mm (±10mm)
Auxiliary Workbench	✓	/	/	/
Keyboard Case	✓	✓	/	✓
Mouse Hanger	✓	/	/	✓
Storage Box	✓	✓	/	✓
Column Material	Aluminum Alloy"	Aluminum Alloy	Aluminum Alloy	Aluminum Alloy
Power System	/	/	Pneumatic Lifting	Pneumatic Lifting
Base	Cast Aluminum with Spraying Coating	Cast Aluminum with Spraying Coating	Cast Aluminum with Spraying Coating	Cast Aluminum with Spraying Coating
Castor	4×4-inch double-sided castors, 2 of which are with brakes	4×4-inch double-sided castors, 2 of which are with brakes	4×4-inch double-sided castors, 2 of which are with brakes	4×4-inch double-sided castors, 2 of which are with brakes
Battery	Optional	Optional	/	Optional

Compatible with All Hospital IT Systems

The workstation cart is fully compatible with HIS, EHR, and PACS systems, supporting seamless digital workflows across departments. It accommodates all leading computer brands and form factors, allowing hospitals to deploy existing IT equipment without additional limitations.

This ensures smooth system integration, faster implementation, and long-term scalability as hospital IT infrastructure evolves.

By combining mobility, ergonomics, modularity, infection control, and IT compatibility, our mobile medical workstation carts are designed to truly improve clinical efficiency and caregiver experience.



SKR-R01

FEATURES AND HIGHLIGHTS



Material

▪ Tabletop:	ABS material
▪ Cart Handle:	One-piece fully cast aluminum, with surface spraying treatment
▪ Keyboard Tray:	ABS material
▪ Drawers:	ABS drawer panels, aluminum alloy drawers
▪ Cart Side Panels and Back Panel:	ABS material
▪ Base:	Q235 material with surface spraying treatment

Technical configuration

☒ Storage Boxes	3 pcs
☒ Mouse Hanger	1 pc
☒ Trash Bins	2 pcs
☒ 4" Double-Sided Silent Castors	4 pcs

Monitor Stand

Detachable Monitor Stand:
Mounting holes (75mm*75mm,
100mm*100mm); constant-force
lifting.




Table Top

ABS medical-grade plastic
tabletop, which is easy to clean
and corrosion-resistant. It has
no surface spraying, so it will not
affect medicines. Grooves are
arranged around the tabletop to
prevent leakage of medical liquids,
and it is also equipped with a
tabletop barrier to prevent items
from falling.



Keyboard Drawer

Concealed keyboard drawer, which
can be stored inside the tabletop
when not in use.

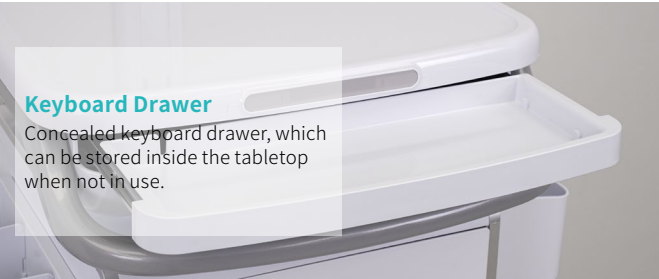


Drawer

4 layers of drawers with a central
control lock, all equipped with self-
closing rails. Internal dimensions
of the first drawer: Length 300
× Width 380 × Height 75 ± 5
mm; Internal dimensions of the
remaining three drawers: Length
300 × Width 380 × Height 165 ±
5 mm.



Handle

Circular anti-collision handles:
These are sturdy metal handles
with rounded corners on all sides,
located at the outermost perimeter
of the entire cart. They not only
prevent injury to other personnel in
case of collisions but also protect
the cart from damage.

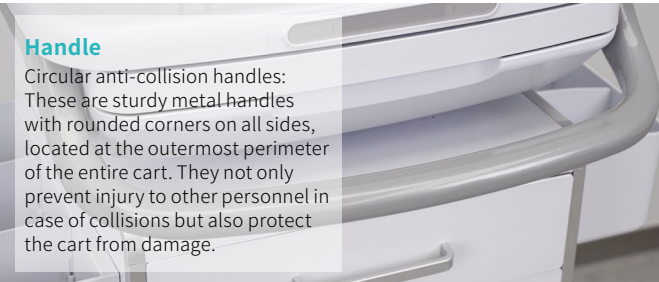

Castors

4-inch double-wheel medical-
grade silent castors, 2 of which are
equipped with brakes.



3 Storage Boxes



1 Mouse Hanger



2 Trash Bins



SKR-R02

FEATURES AND HIGHLIGHTS



Material

▪ Tabletop:	Made of ABS material
▪ Cart Handle:	One-piece fully cast aluminum with surface spraying
▪ Drawers:	ABS drawer panels, aluminum alloy drawer bodies
▪ Columns:	Made of aluminum alloy material
▪ Base:	One-piece fully cast aluminum with surface spraying

Technical configuration

☒ Storage Boxes	2 pcs
☒ Trash Bins	2 pcs
☒ 4" Double-Sided MuteCastors	4 pcs


Monitor Stand

Detachable Monitor Stand:
Mounting holes (75mm*75mm,
100mm*100mm); constant-force
lifting.



Table Top

ABS medical-grade plastic
tabletop, which is easy to clean
and corrosion-resistant. It has
no surface spraying, so it will not
affect medicines. Grooves are
arranged around the tabletop to
prevent leakage of medical liquids,
and it is also equipped with a
tabletop barrier to prevent items
from falling.

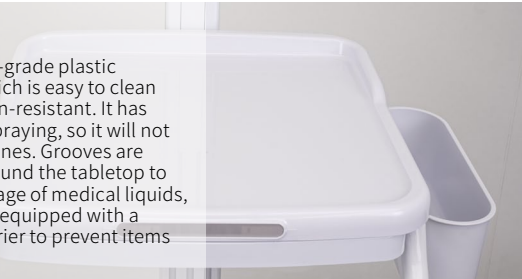

Keyboard Drawer

Concealed keyboard drawer, which
can be stored inside the tabletop
when not in use.



Drawer

Drawers: ABS drawer panels and
aluminum alloy drawer bodies.
Equipped with 2 layers of drawers, a
central control lock, and all drawers
are fitted with self-closing rails.
Internal Dimensions of Drawers:
L300mm*W285mm*H75mm±5mm;




Handle

Circular anti-collision handles:
These are sturdy metal handles
with rounded corners on all sides,
located at the outermost perimeter
of the entire cart. They not only
prevent injury to other personnel in
case of collisions but also protect
the cart from damage.




Column

Material: Aluminum alloy. The
column is equipped with a
professionally calibrated gas
spring, enabling effortless lifting
and lowering. It can be quickly
adjusted to the appropriate height
to accommodate different sitting
and standing postures.

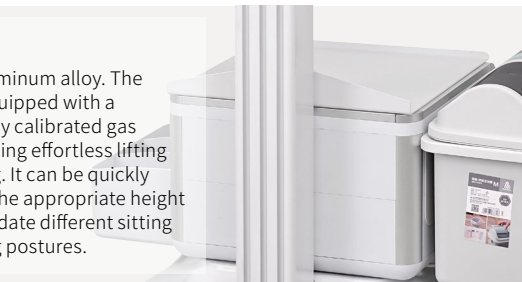


Base

One-piece fully cast aluminum with
surface spraying ensures the center
of gravity and levelness of the
entire cart.

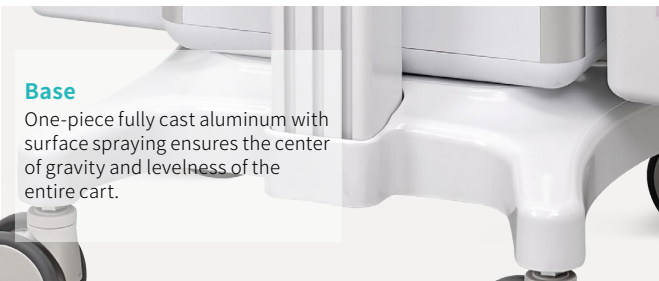

Castors

4 pieces of 4-inch double-sided
mute castors, 2 of which are with
brakes.

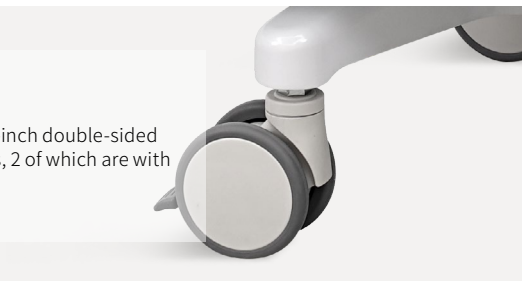

2 Trash Bins



2 Storage Boxes



SKR-R04

FEATURES AND HIGHLIGHTS



Material

- Tabletop: Q235 steel
- Columns: Aluminum alloy material
- Base: One-piece fully cast aluminum with surface spraying

Technical configuration

- ☑ 4" Double-Sided Mute Castors 4 pcs

Monitor Stand

Detachable Monitor Stand:
Mounting holes (75mm*75mm,
100mm*100mm); constant-force
lifting.




Table Top

Material: Q235 steel (with surface
spraying). It is easy to clean
and corrosion-resistant, and is
equipped with stop bars on the left
and right sides to prevent items
from falling.



Lower Tray

The lower tray can hold a printer,
allowing for the easy printing of
patients' medical records and other
materials while conducting mobile
information-based ward rounds.

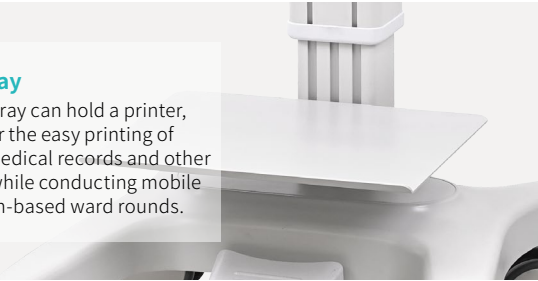

Column

Material: Aluminum alloy. The
column is equipped with a
professionally calibrated gas spring,
enabling effortless lifting and
lowering. It can be quickly adjusted
to the appropriate height to
accommodate different sitting and
standing postures.



Base

One-piece fully cast aluminum with
surface spraying ensures the center
of gravity and levelness of the
entire cart.



Castors

4 pieces of 4-inch double-sided
mute casters, 2 of which are with
brakes.



SKR-R07

FEATURES AND HIGHLIGHTS



Monitor Stand

Detachable Monitor Stand:
Mounting holes (75mm*75mm,
100mm*100mm); constant-force
lifting.




Table Top

ABS medical-grade plastic
tabletop, which is easy to clean
and corrosion-resistant. It has
no surface spraying, so it will not
affect medicines. Grooves are
arranged around the tabletop to
prevent leakage of medical liquids,
and it is also equipped with a
tabletop barrier to prevent items
from falling.



Handle

Circular bumper handles: These are
sturdy metal handles with rounded
corners on all sides, located at
the outermost perimeter of the
entire cart. They not only prevent
injury to other personnel in case of
collisions but also protect the cart
from damage.

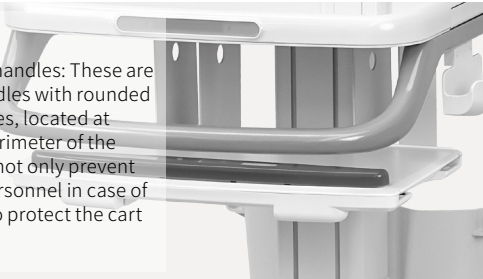


Column

Material: Aluminum alloy. The
column is equipped with a
professionally calibrated gas spring,
enabling effortless lifting and
lowering. It can be quickly adjusted
to the appropriate height to
accommodate different sitting and
standing postures.



Base

One-piece fully cast aluminum with
surface spraying ensures the center
of gravity and levelness of the
entire cart.

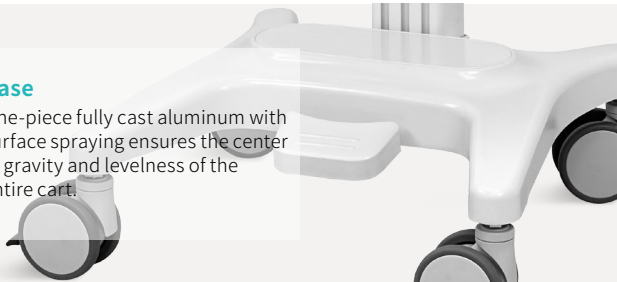

Castors

4 pieces of 4-inch double-sided
mute casters, 2 of which are with
brakes.



One Mouse Hanger



One Storage Box



Material

▪ Tabletop	Made of ABS material
▪ Cart Handle	One-piece fully cast aluminum with surface coating
▪ Columns	Made of aluminum alloy material
▪ Base	One-piece fully cast aluminum with surface coating

Technical configuration

☒ 4" Double-Sided Mute Castors	4 pcs
--	-------

SKR-R05

FEATURES AND HIGHLIGHTS

Technical parameters


- Mainly used for placing tablets;
- Overall Dimensions: L550mm*W550mm*H820-1120mm (±10 mm);
- Tabletop: Carbon steel material, Tabletop Dimensions: L250mm*W300mm (±10 mm);
- Columns: Made of aluminum alloy;
- Lifting Method: Manual lifting, lifting stroke of 300 mm (±10 mm);
- Base: One-piece fully cast aluminum material, with surface painting treatment;
- Casters: 5 pieces of 3-inch double-sided silent casters, 3 of which are equipped with braking function.

Technical configuration

☒ 3" Double-Sided Silent Castors 5 pcs



Table Top
Carbon steel material, Tabletop
Dimensions: Length 250 × Width
300 mm (±10 mm)

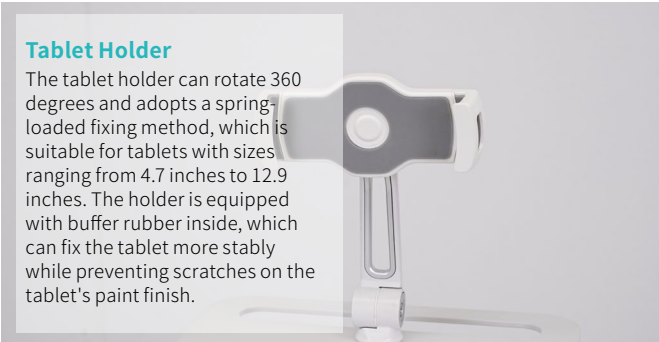


Tablet Holder
The tablet holder can rotate 360
degrees and adopts a spring-
loaded fixing method, which is
suitable for tablets with sizes
ranging from 4.7 inches to 12.9
inches. The holder is equipped
with buffer rubber inside, which
can fix the tablet more stably
while preventing scratches on the
tablet's paint finish.



Column
Aviation aluminum material with
surface spraying treatment;
The column is equipped with a
manual lifting function, allowing
easy upward and downward
movement. It can be quickly
adjusted to the appropriate height
and adopts a knob locking method,
which is suitable for different sitting
and standing postures. The lifting
stroke is 300mm ± 10mm.




Base
One-piece fully cast aluminum
ensures the center of gravity and
levelness of the entire unit, with a
painted surface.

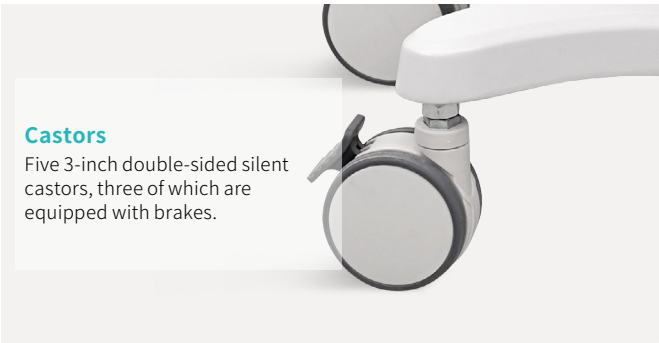

Castors
Five 3-inch double-sided silent
castors, three of which are
equipped with brakes.



CONSTANT INNOVATION

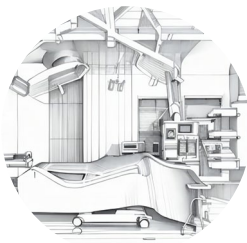


Design



Research



Manufacturing

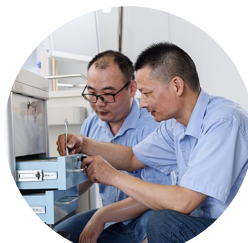

Maintenance



Service



WORLDWIDE BRANCHES

Portugal Showroom

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Brazil Showroom

Rua Lauro Muller 116, Sala
1503 - Botafogo, Rio de
Janeiro - RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

Indian Showroom

1460 M.I.E, Bahadurgarh,
Haryana-124507, INDIA
TEL: +91 96677 70747
E-mail: ind@saikangmedical.
com

China Headquarter

No.35 Lehong Road,
Zhangjiagang City,
Jiangsu Province, China.
TEL: +86 18021231901
E-mail: export@
saikangmedical.com

South East Asia Branch

**SAIKANG Medical
(Thailand) Co., Ltd.**
Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

**Healthward Technology,
INC**
17800 castleton st ste 180,
Los Angeles
TEL: +1 (380) 254 4913

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue